

Your Happy Life
Zen Beam

DESCRIPTION / WARRANTY



MYUNG SHIN MEDICAL CO.,LTD

<http://www.msmedi.com>



Warranty

brand name	ZENBEAM	The term of guarantee
Model	MS-860 CE	one year
Item number		
Day of purchase		

Name of the company	Address	Myung Shin Medical co.,ltd.
		TEL : +82-62-971-8056
Selling agency	35. Cheomdan venture so-ro 62beon-gil, Buk-gu Gwangju 500-460 Korea TEL. +82-62-971-8056 http://www.msmedi.com	

- ※ The Warranty should be kept well as it cannot be reissued.
- ※ Please contact the manufacturer or the selling agency for further information on the warranty.
- ※ This is the topnotch product made by the total quality control and the strict Inspection. In case of a breakdown in normal operating of the product, we will repair it for free under warranty.
- ※ When a breakdown arises from user's carelessness or fixing and remodeling of the product on user's own. You will cover After-sale service expense under warranty.


CONTENTS

1. ZENBEAM ?	4
2. Names of each part of 'Zen-Beam ...	5
3. Names of each part of the internal control	6
4. Manual(Including operating orders) ..	7
5. Directions for use	8
6. How to deal with a breakdown	9
7. Size of the product	10
8. Warranty	attached paper




1 What is 'Zen-Beam' ?

Zen-Beam, which is a thermal equipment conveying Dynamic Wave Energy (Far Infrared Ray) excellent in absorptive power to the body. The product excels in concentrating light by emitting highly efficient Far Infrared Ray in amplifying with heating special Carbon Fiber.

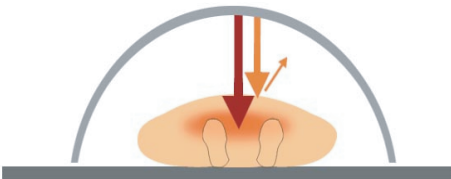


[Principle of Hyper]

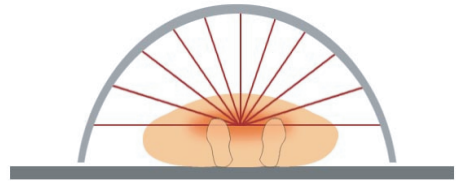


[Principle of concentration]

■ Sectional View of Panel Heater

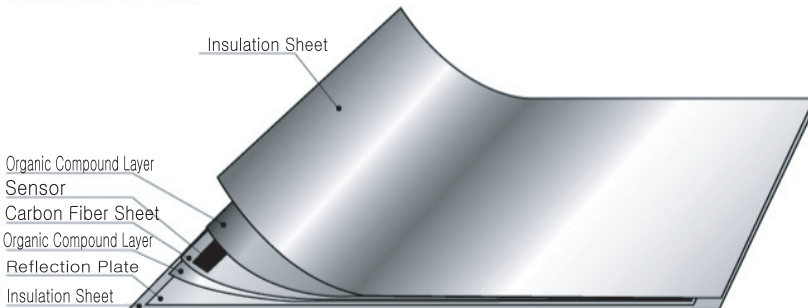




2 Names of each part of 'Zen-Beam'

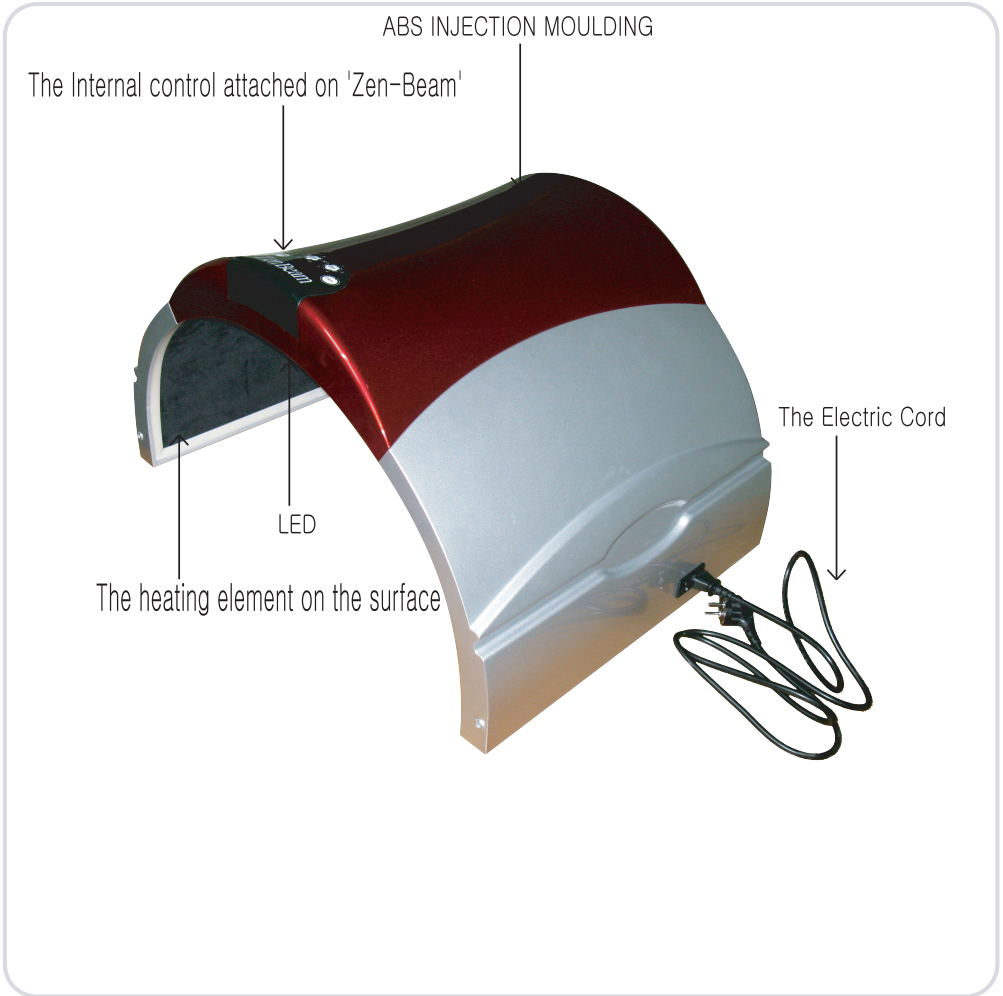



3 Name of each part of the internal control

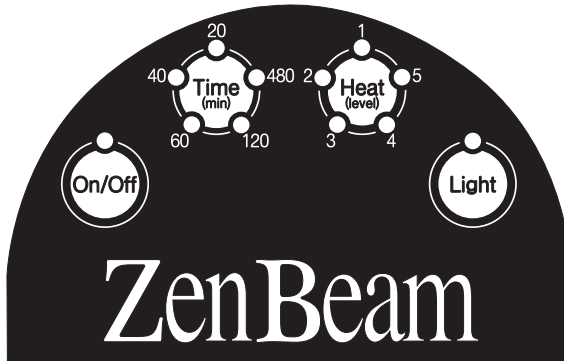

1. ON / OFF : POWER ON / OFF button
2. Time : time control button
3. Heat : heat control button
4. Light : LED on / off button



4 Manual (including operating orders)

- 1) Connect the Electric cable to the outlet
- 2) Button 'ON/OFF'
[First setting Mode : Time – 8hrs(on end), Heat – 5 level]
- 3) Set the proper time pressing Time Button
[Time Control Function : 20min, 40min, 60min, 120min, 8hrs(on end).
Everytime pressing Time Button, setting time is changed.]
- 4) Set the proper temperature pressing Heat Button.
[Heat Control Function : 1~5 level]
- 5) Set the lamp pressing light Button





5 Directions in using

- 1) Do not impact on the product.
Avoid using it at the place there is a lot of moisture.
- 2) Make sure of the voltage for use.
- 3) When using do not put other things on the arch type of the heating plate generating heat.
- 4) Do not take it apart or remodel.
- 5) Do not use it instead of a heater.
- 6) Set the high Heat button because each user has a different constitution.
Watch out for a burn while using the product for a long time.
- 7) Put the cord off the outlet while not using.
- 8) Do not contact with the heating surface.
- 9) Do not use it with the heating surface upside down.
- 10) Please following doctors prescription to use of this product for insensible peripheral nerved diabetic patient, weak skinned women (pregnant), olds and child.
- 11) Please do not use of this product for patients who have an acute disease, heartfailure or set of artificial heart person.
- 12) Please stop to use of this product and following doctors prescription if user get burned during using.
- 13) If there any defects in using, contact the manufacturer of the selling agency.



6 How to deal with breakdown

● When Temperature indicator Lamp Not ON

- Make sure the electric cord and the outlet are connected well.
- Make sure electric supply in the outlet is on.
- Turn on the SWITCH once again.

● When Temperature indicator Lamp Not ON

- After making sure Heat indicator Lamp is ON.
- Turn on the SWITCH once again when Heat indicator Lamp is OFF.

● When Temperature of the Heating Surface Not going up

- Make sure Heat indicator Lamp is ON
- Make sure Time indicator Lamp is ON

※ When STILL NOT ON following making sure of the above items,
Apply to the manufacturer or the selling agency.





7 Size of the product

brand name	ZENBEAM
Model	MS-860 CE
Electricity supply for input	220V~240V
Frequency	50~60Hz
Maximum Power consumption	350W
Size (mm)	651mm(wide) * 480mm(deep) * 380mm(high)
Weight	About 6kg



© The meaning of Symbol Mark ©

The four elements of life are water, air, light and salt. The human body can even purify them by the body's absorbing Far-Infrared ray(light of life) unpolluted. It means that our mark forms the process of absorbing Far-Infrared ray emitted in space into the body and strongly leaps for a global company.





(주)명신메디칼
MYUNG SHIN MEDICAL CO.,LTD

Manufacturer

35, Cheomdan venture so-ro 62beon-gil, Buk-gu
Gwangju 500-460 Korea
TEL.+82-62-971-8056 <http://www.msmedi.com>