


HL7 CONFORMANCE STATEMENT

cv42®

Version 6.5.X
April 30, 2026

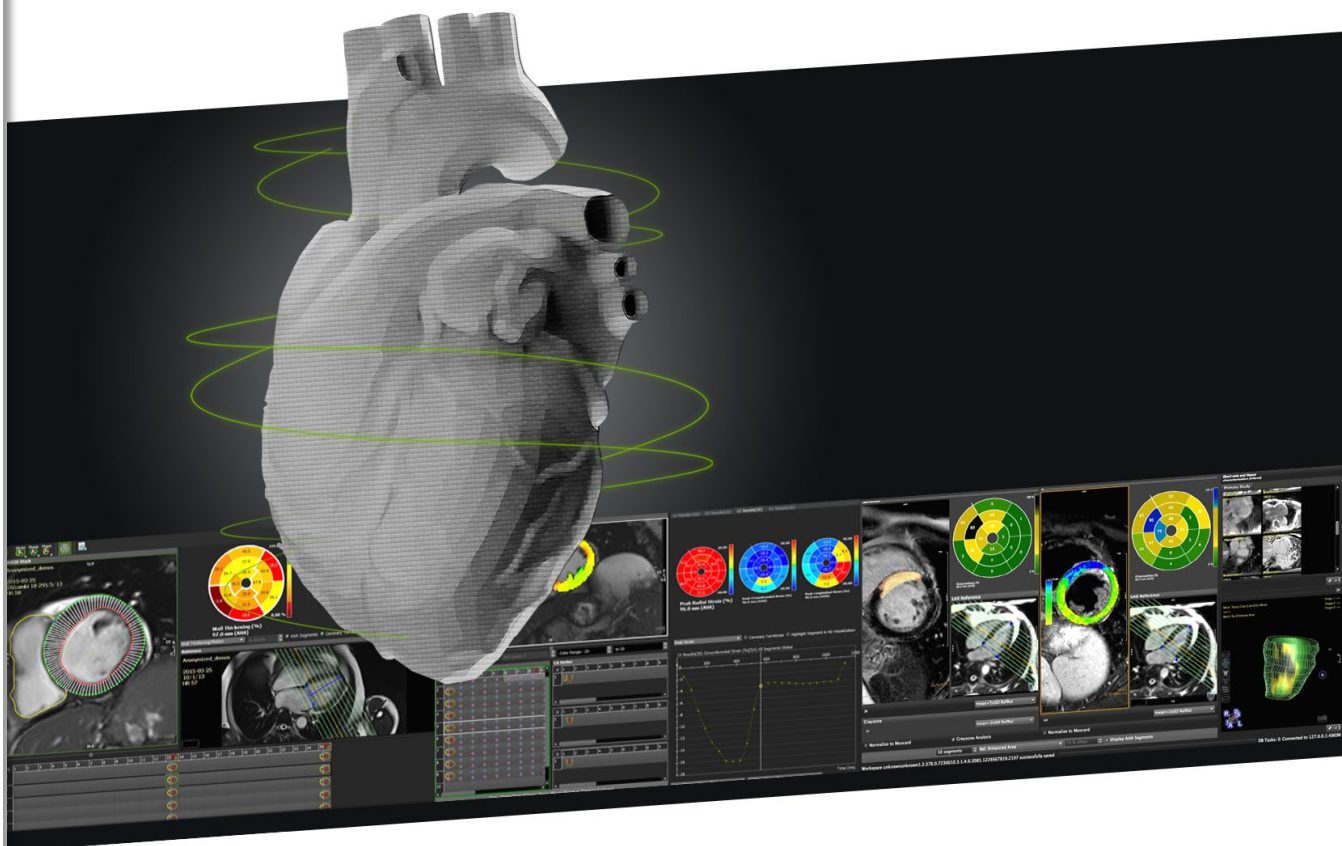

At the **Heart** of
IMAGING

1 CONFORMANCE STATEMENT OVERVIEW

cvi42 reporting is a vendor-independent workstation solution for comprehensive reporting of Cardiovascular MR and CT cases. It is designed to be used in clinical settings as well as for experimental and clinical research.

cvi42 reporting supports HL7 v2.3 ORM (received) and ORU (sent) message, using Mirth Connect as a proxy server to adapt and transform HL7 messages as they are exchanged with an upstream HIS, RIS or HL7 Engine.

2 TABLE OF CONTENTS

1	CONFORMANCE STATEMENT OVERVIEW	2
2	TABLE OF CONTENTS.....	3
3	INTRODUCTION	4
3.1	Revision History	4
3.2	Audience	5
3.3	Remarks.....	5
3.4	Use Case Model	5
3.5	Dynamic Interaction Model	6
3.6	Dynamic Definitions	7
3.6.1	Optionality Codes.....	7
3.6.2	ORM (Event O01)	7
3.6.3	ORU (Event R01).....	7
3.7	Static Definitions – Message Level.....	7
3.7.1	ORM^O01	7
3.7.2	ORU^R01.....	7
3.8	Static Definitions – Segment Level.....	8
3.8.1	MSH	8
3.8.2	PID	9
3.8.3	PV1	9
3.8.4	PV2	9
3.8.5	ORC	10
3.8.6	OBR	11
3.8.7	OBX	13
3.8.8	OBX-3 Observation Identifiers for Discrete Data	15

3 INTRODUCTION

This HL7 Conformance Statement specifies the behavior and functionality of the cvi42 reporting system, with regard to supported HL7 Messages and Message Segments. cvi42 reporting is an intranet web application for comprehensive evaluation and reporting of Cardiovascular MR and CT data.

Contact address:

CIRCLE Cardiovascular Imaging Inc
Suite 1800, 707 – 8th Ave SW
Calgary, Alberta
Canada T2P 1H5
<http://www.circlecvi.com/>

3.1 Revision History

Document Version	Date of Issue	Author	Description
Version report42 3.0.4	2015-10-06	Glen van de Mosselaer	Initial release of this document
Version report42 3.1.1	2016-02-22	Glen van de Mosselaer	Revision for version 3.1.x functionality
Version report42 3.2.0	2016-08-09	Glen van de Mosselaer	Revision for version 3.2.x functionality
Version report42 3.3.0	2017-02-01	Glen van de Mosselaer	Revision for version 3.3.x functionality
Version report42 3.5.0	2018-03-28	Angad Thandi	Revision for version 3.5.x functionality
Version report42 3.6.0	2018-11-28	Billy Wu	Revision for version 3.6.x functionality
Version report42 3.7.0	2019-07-25	Billy Wu	Version update
Version report42 3.8.0	2020-04-01	Angad Thandi	Revision for version 3.8.x functionality
Version report42 3.9.0	2020-07-22	Angad Thandi	Revision for version 3.9.x functionality
Version report42 3.9.5	2021-03-08	Bill Luu	Revision for version 3.9.5 functionality
Version report42 3.9.5	2021-04-08	Bill Luu	Revision for correct Alternate Identifier OBX-3.4 values
Version report42 3.9.6	2021-04-26	Andrew Schneider	Revision for Myocardial Strain OBX-3.2 values
Version report42 3.10.0	2022-03-08	Bill Luu	Revision for version 3.10.x functionality
Version report 3.10.1	2022-05-02	Andrew Schneider	Revision for version 3.10.1 functionality
Version report 3.10.1	2022-06-16	Bill Luu	Version update
Version report 5.15.0	2022-10-28	Angad Thandi	Revision for version 5.15.0 functionality
Version report 5.16.0	2022-12-13	Nishanth Gandhi	Update to section 3.5 Dynamic Interaction Model

Version report 5.17.0	2023-03-27	Angad Thandi	Revision for version 5.17.0 functionality
Version report 6.0.0	2023-10-05	Angad Thandi	Revision for version 6.0.0 functionality
Version report 6.1.0	2024-04-16	Angad Thandi	Revision for version 6.1.0 functionality
Version cvi42 6.1.0	2024-04-16	Derek Hann	Updating document to reflect that reporting is now a part of cvi42
Version cvi42 6.2.0	2024-12-05	Angad Thandi	Revision for version 6.2.0 functionality
Version cvi42 6.3.0	2025-05-23	Angad Thandi	Revision for version 6.3.0 functionality
Version cvi42 6.3.1	2025-07-22	Angad Thandi	Revision for version 6.3.1 functionality
Version cvi42 6.4.0	2025-11-03	Angad Thandi	Revision for version 6.4.0 functionality
Version cvi42 6.5.0	2026-04-30	Angad Thandi	Revision for version 6.5.0 functionality

3.2 Audience

This document is written for the people that need to understand how **cvi42** reporting will integrate into their healthcare facility. This includes both those responsible for overall workflow network policy and architecture, as well as integrators who need to have a detailed understanding of the HL7 features of the product. Integrators are expected to fully understand all the HL7 terminology, how the tables in this document relate to the product's functionality, and how that functionality integrates with other devices that support compatible HL7 features.

3.3 Remarks

The scope of this HL7 Conformance Statement is to facilitate integration between **cvi42** reporting and other HL7 products. The Conformance Statement should be read and understood in conjunction with the HL7 v2.3 Standard. HL7 by itself does not guarantee interoperability. The Conformance Statement does, however, facilitate a first-level comparison for interoperability between different applications supporting compatible HL7 functionality.

This Conformance Statement is not supposed to replace validation with other HL7 equipment to ensure proper exchange of intended information. In fact, the user should be aware of the following important issues:


- The comparison of different Conformance Statements is just the first step towards assessing interconnectivity and interoperability between the product and other HL7 conformant equipment.
- Test procedures should be defined and executed to validate the required level of interoperability with specific compatible HL7 equipment, as established by the healthcare facility.
- Whole the message content (XML) is saved to the database in its raw-format in tblHI7ImagingOrder

3.4 Use Case Model

cvi42 reporting can receive ORM (Order Message) from a HIS, RIS, or general HL7 Engine indicating that a patient was scheduled for an MR or CT Scan and the order is complete. These messages should be received prior to actual image or measurement data being transferred from **cvi42** or other analysis software.

When a report is changed to a Finalized or Addendum state, an ORU (Observation Message) is returned to an upstream HIS, RIS, or HL7 Engine, presumably but not necessarily the same system that originated the transaction with an ORM.

3.5 Dynamic Interaction Model

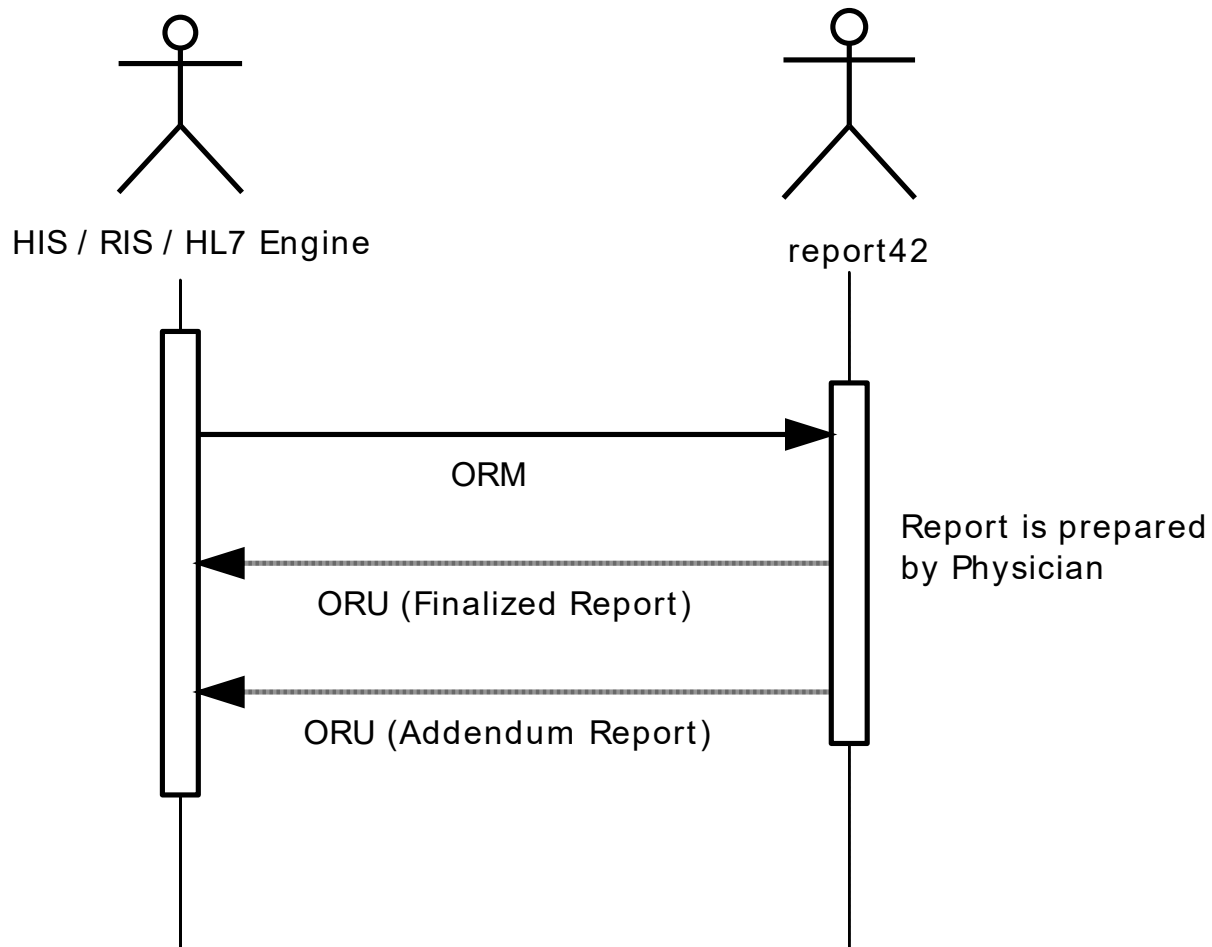

ORM messages originate with a HIS, RIS or HL7 Engine upstream from the cvi42 reporting server. Received ORM messages are stored in cvi42 reporting's database and matched later to incoming DICOM-sourced data.

ORU messages are later produced when a physician creates a Finalized or Addendum report and are returned to an HL7 node, presumably the one that created the ORM message. The nodes and the messages are matched in the following order:

- StudyUID – the unique id the study has in the reporting database
- Accession Number

Typically report data will be included either by embedding text or encoded data, or by including reference pointers to report data at a local or network location.

3.6 Dynamic Definitions

3.6.1 Optionality Codes

R: Required
 RE: Required but may be empty
 O: Optional
 C: Conditional
 X: Not supported
 B: Backward Compatibility
 W: Withdrawn

3.6.2 ORM (Event O01)

MSH PID PV1 [PV2] ORC [OBR] [{ OBX }]

3.6.3 ORU (Event R01)

MSH PID PV1 [PV2] ORC OBR [{ OBX }]

3.7 Static Definitions – Message Level

3.7.1 ORM^O01

Segment	ORM Message	Usage	Cardinality	Comment
MSH	Message Header	R	[1..1]	
PID	Patient Identification	R	[1..1]	
PV1	Patient Visit	RE	[1..1]	
[PV2]	Patient Visit – Additional	O	[0..1]	
ORC	Common Order	R	[1..1]	Required in order to filter messages by Order Status
[OBR]	Observation	O	[0..1]	
[{ OBX }]	Observation/Result	O	[0..*]	

3.7.2 ORU^R01

Segment	ORM Message	Usage	Cardinality	Comment
MSH	Message Header	R	[1..1]	
PID	Patient Identification	R	[1..1]	
PV1	Patient Visit	RE	[1..1]	
[PV2]	Patient Visit – Additional	O	[0..1]	If present in ORM, will be returned in ORU
ORC	Common Order	R	[1..1]	
OBR	Observation	O	[1..1]	
[{ OBX }]	Observation/Result	O	[0..*]	OBX present in ORM will not be present in ORU. Only specified report data will be returned in OBX segment(s)

3.8 Static Definitions – Segment Level

3.8.1 MSH

Sequence	Name	Usage	Cardinality	Comment
1	Field Separator	R	[1..1]	
2	Encoding Characters	R	[1..1]	
3	Sending Application	R	[1..1]	Sender in ORM becomes Receiver in ORU. Interpreted
4	Sending Facility	R	[1..1]	Sender in ORM becomes Receiver in ORU. Interpreted
5	Receiving Application	O	[0..1]	
6	Receiving Facility	O	[0..1]	
7	Date/Time Of Message	R	[0..1]	Interpreted
8	Security	O	[0..1]	
9	Message Type	R	[1..1]	
10	Message Control ID	R	[1..1]	
11	Processing ID	R	[1..1]	Interpreted
12	Version ID	R	[1..1]	Interpreted
13	Sequence Number	O	[0..1]	
14	Continuation Pointer	O	[0..1]	
15	Accept Acknowledgement Type	O	[0..1]	
16	Application Acknowledgement Type	O	[0..1]	
17	Country Code	O	[0..1]	
18	Character Set	O	[0..*]	
19	Principal Language Of Message	O	[0..1]	
20	Alternate Character Set Handling Scheme	O	[0..1]	

3.8.2 PID

Sequence	Name	Usage	Cardinality	Comment
1	Set ID – PID	O	[0..1]	
2	Patient ID	B	[0..1]	Interpreted
3	Patient Identifier List	R	[1..*]	Interpreted
4	Alternate Patient ID – PID	B	[1..*]	
5	Patient's Name	R	[1..*]	Interpreted
6	Mother's Maiden Name	O	[0..*]	
7	Date/Time of Birth	O	[0..1]	Interpreted
8	Sex	O	[0..1]	Interpreted
9	Patient Alias	O	[0..*]	
10	Race	O	[0..*]	
11	Patient Address	O	[0..*]	
12	Country Code	B	[0..1]	
13	Phone Number – Home	O	[0..*]	
14	Phone Number – Business	O	[0..*]	
15	Primary Language	O	[0..1]	
16	Marital Status	O	[0..1]	
17	Religion	O	[0..1]	
18	Patient Account Number	O	[0..1]	
19	SSN Number – Patient	B	[0..1]	
20	Driver's License Number – Patient	O	[0..1]	
21	Mother's Identifier	O	[0..*]	
22	Ethnic Group	O	[0..*]	
23	Birth Place	O	[0..1]	
24	Multiple Birth Indicator	O	[0..1]	
25	Birth Order	O	[0..1]	
26	Citizenship	O	[0..*]	
27	Veterans Military Status	O	[0..1]	
28	Nationality	O	[0..1]	
29	Patient Death Date and Time	O	[0..1]	
30	Patient Death Indicator	O	[0..1]	

3.8.3 PV1

PV1 segments are ignored by cvi42 reporting, however the segment will be stored as-is from an ORM and returned in the ORU.

3.8.4 PV2

PV2 segments are ignored by cvi42 reporting, however the segment will be stored as-is from an ORM and returned in the ORU.

3.8.5 ORC

Sequence	Name	Usage	Cardinality	Comment
1	Order Control	R	[1..1]	
2	Placer Order Number	O	[0..1]	
3	Filler Order Number	O	[0..1]	
4	Placer Group Number	O	[0..1]	
5	Order Status	O	[0..1]	
6	Response Flag	O	[0..1]	
7	Quantity/Timing	O	[0..1]	
8	Parent	O	[0..1]	
9	Date/Time of Transaction	O	[0..1]	
10	Entered By	O	[0..*]	
11	Verified By	O	[0..*]	
12	Ordering Provider	O	[0..*]	
13	Enterer's Location	O	[0..1]	
14	Call Back Phone Number	O	[0..2]	
15	Order Effective Date/Time	O	[0..1]	
16	Order Control Code Reason	O	[0..1]	
17	Entering Organization	O	[0..1]	
18	Entering Device	O	[0..1]	
19	Action By	O	[0..*]	
20	Advanced Beneficiary Notice Code	O	[0..1]	
21	Ordering Facility Name	O	[0..*]	
22	Ordering Facility Address	O	[0..*]	
23	Ordering Facility Phone Number	O	[0..*]	
24	Ordering Provider Address	O	[0..*]	

3.8.6 OBR

Sequence	Name	Usage	Cardinality	Comment
1	Set ID – OBR	O	[0..1]	
2	Placer Order Number	O	[0..1]	
3	Filler Order Number	O	[0..1]	
4	Universal Service ID	R	[1..1]	
5	Priority	O	[0..1]	
6	Requested Date/Time	O	[0..1]	
7	Observation Date/Time	O	[0..1]	
8	Observation End Date/Time	O	[0..1]	
9	Collection Volume	O	[0..1]	
10	Collector Identifier	O	[0..*]	
11	Specimen Action Code	O	[0..1]	
12	Danger Code	O	[0..1]	
13	Relevant Clinical Info.	O	[0..1]	
14	Specimen Received Date/Time	O	[0..1]	
15	Specimen Source	O	[0..1]	
16	Ordering Provider	O	[0..*]	
17	Order Callback Phone Number	O	[0..2]	
18	Placer Field 1	O	[0..1]	
19	Placer Field 2	O	[0..1]	
20	Filler Field 1	O	[0..1]	For ORU will contain the DICOM StudyInstanceUID
21	Filler Field 2	O	[0..1]	
22	Results Rpt/Status Chng – Date/Time	O	[0..1]	For ORU will contain the date/time of the message
23	Charge to Practice	O	[0..1]	
24	Diagnostic Serv Sect ID	O	[0..1]	
25	Result Status	O	[0..1]	In ORU message may be 'F' for Finalized or 'C' for Correction (Addendum) reports
26	Parent Status	O	[0..1]	
27	Quantity/Timing	O	[0..*]	
28	Result Copies To	O	[0..5]	
29	Parent	O	[0..1]	
30	Transportation Mode	O	[0..1]	
31	Reason for Study	O	[0..*]	
32	Principal Result Interpreter	O	[0..1]	In ORU message set to the reporting physician
33	Assistant Result Interpreter	O	[0..*]	
34	Technician	O	[0..*]	
35	Transcriptionist	O	[0..*]	
36	Scheduled Date/Time	O	[0..1]	
37	Number of Sample Containers	O	[0..1]	
38	Transport Logistics of Collected Sample	O	[0..*]	
39	Collector's Comment	O	[0..*]	
40	Transport Arrangement Responsibility	O	[0..1]	
41	Transport Arranged	O	[0..1]	
42	Escort Required	O	[0..1]	

43	Planned Patient Transport Comment	O	[0..*]	
44	Procedure Code	O	[0..1]	
45	Procedure Code Modifier	O	[0..*]	

3.8.7 OBX

Sequence	Name	Usage	Cardinality	Comment
1	Set ID – OBX	O	[0..1]	
2	Value Type	O	[0..1]	ORU messages use TX for Text report data, ED for Base64-encoded PDF data, and RP for references to report files stored on a filesystem or a DICOM SeriesInstanceUID for reports uploaded to PACS. For Discrete Data, this may be TX for blocks of Text data, ST for a String (short alphanumeric text), or NM for Numeric data.
3	Observation Identifier	R	[1..1]	For embedded formatted reports or formatted report reference pointers, subfield 1 contains a LOINC code ID, subfield 2 contains a LOINC Long Common Name, and subfield 3 contains LN signifying that this observation is from a Heart CT (58744-4^Heart CT^LN) or Heart MRI (24748-6^Heart MRI^LN). For discrete data points, subfield 1 contains a short data identifier, subfield 2 contains more descriptive identifier text, and subfield 3 contains the name of the cvi42 reporting coding system, 99R42. For discrete data points with relevant LOINC codes, the LOINC code will be represented in fields 4, 5, and 6. Refer to the table of OBX-3 Observation Identifiers for Discrete Data for a complete list of discrete values and codes.
4	Observation Sub-ID	O	[0..1]	For Discrete Data, may indicate that a data point belongs to a particular region of interest or slices (OBX-2 Value Types ST or NM), or may indicate individual segments of a Text value (OBX-2 Value Type TX).
5	Observation Value	O	[0..*]	For embedded Text, there are no subfields and OBX-5 contains the embedded content.

				For Embedded PDF, subfield 1 contains "Adobe Acrobat", subfield 2 contains "PDF", subfield 4 contains "Base64" indicating that this data is encapsulated using Base64 encoding, and subfield 5 contains the actual data. For Reference Pointers, subfield 1 contains the path or URL to the document, subfield 3 contains the name of the data type ('text' for text, 'application' for PDF, 'image' for DICOM) and subfield 4 contains the MIME type ('text/plain', 'application/pdf', or 'image/dicom'). For Discrete Data, the field will contain the actual value of the data point.
6	Units	O	[0..1]	For Discrete Data, will contain the units for the data, if any.
7	References Range	O	[0..1]	For Discrete Data, will contain the reference (normal value) range for the data, if any.
8	Abnormal Flags	O	[0..5]	For Discrete Data, will contain a flag indicating if the data is in-range (N), higher than the reference range (H) or lower than the reference range (L), if a reference range is defined.
9	Probability	O	[0..1]	
10	Nature of Abnormal Test	O	[0..*]	
11	Observation Result Status	R	[1..1]	In ORU message may be 'F' for Finalized or 'C' for Correction (Addendum) reports
12	Date Last Obs Normal Values	O	[0..1]	
13	User Defined Access Checks	O	[0..1]	
14	Date/Time of the Observation	O	[0..1]	
15	Producer's ID	O	[0..1]	
16	Responsible Observer	O	[0..*]	
17	Observation Method	O	[0..*]	

3.8.8 OBX-3 Observation Identifiers for Discrete Data

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
Patient / Study					
SiteOfServiceLaboratory	Site Of Service Laboratory Name	99R42			
SiteOfServiceLocation	Site Of Service Location	99R42			
SiteOfServiceTypeOfFacility	Site Of Service Type Of Facility	99R42			
SiteOfServiceAddress	Site Of Service Address	99R42			
SiteOfServicePhone	Site Of Service Phone	99R42			
SiteOfServiceAccreditationEntity	Site Of Service Accreditation Entity	99R42			
SiteOfServiceAccreditationStatus	Site Of Service Accreditation Status	99R42			
PHN	Patient Health Number	99R42			
StudyInstanceUID	DICOM StudyInstanceUID	99R42			
AccessionNumber	DICOM Accession Number	99R42			
PriorNames	Patient Prior Names	99R42			
Ethnicity	Patient Ethnicity	99R42			
InsuranceInformation	Patient Insurance Information	99R42			
Status	Report Status	99R42			
MRN	Medical Record Number	99R42			
CustomPatientID	DICOM Custom Patient ID	99R42			
Gender	Patient Gender	99R42			
Birthdate	Patient Birthdate	99R42			
StudyDate	DICOM Study Date	99R42			
PatientAge	Patient Age	99R42			
StudyName	DICOM Study Name	99R42			
ReferringPhysician	DICOM Referring Physician	99R42			
PhysicianNationalIdentifier	Referring Physician National Identifier	99R42			
BloodPressure	Blood Pressure	99R42			
HeartRate	Heart Rate	99R42			
Height	Height	99R42			
Weight	Weight	99R42			
BMI	Body Mass Index	99R42			
BSA	Body Surface Area	99R42			
BSAMethod	Body Surface Area Formula	99R42			
ImageQuality	Image Quality	99R42			
ContrastAgentAdministrationMethod	Contrast Agent Administration Method	99R42			
ContrastAgent	Contrast Agent	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
ContrastAgentDose	Contrast Agent Dose	99R42			
ScannerManufacturer	Scanner Manufacturer	99R42			
ScannerModel	Scanner Model	99R42			
ScannerSerialNumber	Scanner Serial Number	99R42			
ScannerSoftwarePlatform	Scanner Software Platform	99R42			
InducingStressMethod	Inducing Stress Method	99R42			
StressProtocol	Stress Protocol	99R42			
MagnetType	Magnet Type	99R42			
FieldStrength	Field Strength	99R42			
ScannerType	Scanner Type	99R42			
Gating	Gating	99R42			
mAs	mAs	99R42			
kVp	kVp	99R42			
RadiationReductionStrategy	Radiation Reduction Strategy	99R42			
EstimatedRadiationDose	Estimated Radiation Dose	99R42			
PatientHistory	Patient History	99R42			
ClinicalSituations	Clinical Situations	99R42			
StudyIndications	Study Indications	99R42			
RiskFactors	Risk Factors	99R42			
Medications	Medications	99R42			
Allergies	Allergies	99R42			
Staff	Staff	99R42			
Modality	Modality	99R42			
ProtocolName	Protocol Name	99R42			
IndicationName	Indication Name	99R42			
Techniques	Techniques	99R42			
LV Global					
LV_EDV	LV ED volume	99R42	8821-1	LV ED vol Imaging	LN
LV_EDV_ZValue	LV ED volume Z-Value	99R42			
LV_EDVI	LV ED volume / height	99R42			
LV_EDV_AI	LV ED volume / BSA	99R42			
LV_ESV	LV ES volume	99R42	8823-7	LV ES vol Imaging	LN
LV_ESV_ZValue	LV ES volume Z-Value	99R42			
LV_ESVI	LV ES volume / Height	99R42			
LV_ESV_AI	LV ES volume / BSA	99R42			
LV_SV	LV SV	99R42	8772-6	LV SV MRI	LN
LV_SV_ZValue	LV SV Z-Value	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
LV_SVI	LV SV / height	99R42			
LV_SV_AI	LV SV / BSA	99R42	8792-4	LV SVI MRI	LN
LV_EF	LV EF	99R42	8811-2	LV EF MRI	LN
LV_EF_ZValue	LV EF Z-Value	99R42			
LV_CardiacOutput	LV Cardiac Output	99R42	8738-7	LV Output MRI	LN
LV_CardiacOutput_ZValue	LV Cardiac Output Z-Value	99R42			
LV_CardiacIndex	LV Cardiac Index	99R42	8752-8	LV Cardiac index MRI	LN
LV_Mass	LV Mass	99R42			
LV_Mass_ZValue	LV Mass Z-Value	99R42			
LV_MI	LV Mass / height	99R42			
LV_Mass_AI	LV Mass / BSA	99R42			
LV_ESWS	LV ES Wall Stress	99R42			
LVGlobalMass_Method	LV Global Mass Method	99R42			
LV_MyoMass_Diast	LV Myo Mass Diastole	99R42			
LV_MyoMass_Diast_I	LV Myo Mass Diastole / height	99R42			
LV_MyoMass_Diast_AI	LV Myo Mass Diastole / BSA	99R42			
LV_MyoMass_Syst	LV Myo Mass Systole	99R42			
LV_MyoMass_Syst_I	LV Myo Mass Systole / height	99R42			
LV_MyoMass_Syst_AI	LV Myo Mass Systole / BSA	99R42			
LV_PapMass_Diast	LV Papillary Mass Diastole	99R42			
LV_PapMass_Syst	LV Papillary Mass Systole	99R42			
LV_HeartRate	LV Heart Rate	99R42			
GlobalPeakWallThickness_diast	Global Peak Wall Thickness	99R42			
MAPSE_Inferior	MAPSE Inferior	99R42			
MAPSE_Anterior	MAPSE Anterior	99R42			
MAPSE_Lateral	MAPSE Lateral	99R42			
MAPSE_Septal	MAPSE Septal	99R42			
LV Regional Tissue Characterization					
LVRegTissueChar_Endocardium_Segment_1	LV Regional Tissue Characterization: Endocardium Segment 1	99R42			
LVRegTissueChar_Epicardium_Segment_1	LV Regional Tissue Characterization: Epicardium Segment 1	99R42			
LVRegTissueChar_Endocardium_Segment_2	LV Regional Tissue Characterization: Endocardium Segment 2	99R42			
LVRegTissueChar_Epicardium_Segment_2	LV Regional Tissue Characterization: Epicardium Segment 2	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
LVRegTissueChar_Endocardium_Segment_3	LV Regional Tissue Characterization: Endocardium Segment 3	99R42			
LVRegTissueChar_Epicardium_Segment_3	LV Regional Tissue Characterization: Epicardium Segment 3	99R42			
LVRegTissueChar_Endocardium_Segment_4	LV Regional Tissue Characterization: Endocardium Segment 4	99R42			
LVRegTissueChar_Epicardium_Segment_4	LV Regional Tissue Characterization: Epicardium Segment 4	99R42			
LVRegTissueChar_Endocardium_Segment_5	LV Regional Tissue Characterization: Endocardium Segment 5	99R42			
LVRegTissueChar_Epicardium_Segment_5	LV Regional Tissue Characterization: Epicardium Segment 5	99R42			
LVRegTissueChar_Endocardium_Segment_6	LV Regional Tissue Characterization: Endocardium Segment 6	99R42			
LVRegTissueChar_Epicardium_Segment_6	LV Regional Tissue Characterization: Epicardium Segment 6	99R42			
LVRegTissueChar_Endocardium_Segment_7	LV Regional Tissue Characterization: Endocardium Segment 7	99R42			
LVRegTissueChar_Epicardium_Segment_7	LV Regional Tissue Characterization: Epicardium Segment 7	99R42			
LVRegTissueChar_Endocardium_Segment_8	LV Regional Tissue Characterization: Endocardium Segment 8	99R42			
LVRegTissueChar_Epicardium_Segment_8	LV Regional Tissue Characterization: Epicardium Segment 8	99R42			
LVRegTissueChar_Endocardium_Segment_9	LV Regional Tissue Characterization: Endocardium Segment 9	99R42			
LVRegTissueChar_Epicardium_Segment_9	LV Regional Tissue Characterization: Epicardium Segment 9	99R42			
LVRegTissueChar_Endocardium_Segment_10	LV Regional Tissue Characterization: Endocardium Segment 10	99R42			
LVRegTissueChar_Epicardium_Segment_10	LV Regional Tissue Characterization: Epicardium Segment 10	99R42			
LVRegTissueChar_Endocardium_Segment_11	LV Regional Tissue Characterization: Endocardium Segment 11	99R42			
LVRegTissueChar_Epicardium_Segment_11	LV Regional Tissue Characterization: Epicardium Segment 11	99R42			
LVRegTissueChar_Endocardium_Segment_12	LV Regional Tissue Characterization: Endocardium Segment 12	99R42			
LVRegTissueChar_Epicardium_Segment_12	LV Regional Tissue Characterization: Epicardium Segment 12	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
LVRegTissueChar_Endocardium_Segment_13	LV Regional Tissue Characterization: Endocardium Segment 13	99R42			
LVRegTissueChar_Epicardium_Segment_13	LV Regional Tissue Characterization: Epicardium Segment 13	99R42			
LVRegTissueChar_Endocardium_Segment_14	LV Regional Tissue Characterization: Endocardium Segment 14	99R42			
LVRegTissueChar_Epicardium_Segment_14	LV Regional Tissue Characterization: Epicardium Segment 14	99R42			
LVRegTissueChar_Endocardium_Segment_15	LV Regional Tissue Characterization: Endocardium Segment 15	99R42			
LVRegTissueChar_Epicardium_Segment_15	LV Regional Tissue Characterization: Epicardium Segment 15	99R42			
LVRegTissueChar_Endocardium_Segment_16	LV Regional Tissue Characterization: Endocardium Segment 16	99R42			
LVRegTissueChar_Epicardium_Segment_16	LV Regional Tissue Characterization: Epicardium Segment 16	99R42			
LVRegTissueChar_Epicardium_Segment_17	LV Regional Tissue Characterization: Epicardium Segment 17	99R42			
LVSegFunc_Segment_1	LV Regional Segmental Function: Segment 1	99R42			
LVSegFunc_Segment_2	LV Regional Segmental Function: Segment 2	99R42			
LVSegFunc_Segment_3	LV Regional Segmental Function: Segment 3	99R42			
LVSegFunc_Segment_4	LV Regional Segmental Function: Segment 4	99R42			
LVSegFunc_Segment_5	LV Regional Segmental Function: Segment 5	99R42			
LVSegFunc_Segment_6	LV Regional Segmental Function: Segment 6	99R42			
LVSegFunc_Segment_7	LV Regional Segmental Function: Segment 7	99R42			
LVSegFunc_Segment_8	LV Regional Segmental Function: Segment 8	99R42			
LVSegFunc_Segment_9	LV Regional Segmental Function: Segment 9	99R42			
LVSegFunc_Segment_10	LV Regional Segmental Function: Segment 10	99R42			
LVSegFunc_Segment_11	LV Regional Segmental Function: Segment 11	99R42			
LVSegFunc_Segment_12	LV Regional Segmental Function: Segment 12	99R42			
LVSegFunc_Segment_13	LV Regional Segmental Function: Segment 13	99R42			
LVSegFunc_Segment_14	LV Regional Segmental Function: Segment 14	99R42			
LVSegFunc_Segment_15	LV Regional Segmental Function: Segment 15	99R42			
LVSegFunc_Segment_16	LV Regional Segmental Function: Segment 16	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
LVSegFunc_Segment_17	LV Regional Segmental Function: Segment 17	99R42			
LVSegViability_Segment_1	LV Regional Segmental Viability: Segment 1	99R42			
LVSegViability_Segment_2	LV Regional Segmental Viability: Segment 2	99R42			
LVSegViability_Segment_3	LV Regional Segmental Viability: Segment 3	99R42			
LVSegViability_Segment_4	LV Regional Segmental Viability: Segment 4	99R42			
LVSegViability_Segment_5	LV Regional Segmental Viability: Segment 5	99R42			
LVSegViability_Segment_6	LV Regional Segmental Viability: Segment 6	99R42			
LVSegViability_Segment_7	LV Regional Segmental Viability: Segment 7	99R42			
LVSegViability_Segment_8	LV Regional Segmental Viability: Segment 8	99R42			
LVSegViability_Segment_9	LV Regional Segmental Viability: Segment 9	99R42			
LVSegViability_Segment_10	LV Regional Segmental Viability: Segment 10	99R42			
LVSegViability_Segment_11	LV Regional Segmental Viability: Segment 11	99R42			
LVSegViability_Segment_12	LV Regional Segmental Viability: Segment 12	99R42			
LVSegViability_Segment_13	LV Regional Segmental Viability: Segment 13	99R42			
LVSegViability_Segment_14	LV Regional Segmental Viability: Segment 14	99R42			
LVSegViability_Segment_15	LV Regional Segmental Viability: Segment 15	99R42			
LVSegViability_Segment_16	LV Regional Segmental Viability: Segment 16	99R42			
LVSegViability_Segment_17	LV Regional Segmental Viability: Segment 17	99R42			
RV Global					
RV_EDV	RV ED volume	99R42	8822-9	RV ED vol Imaging	LN
RV_EDV_ZValue	RV ED volume Z-Value	99R42			
RV_EDVI	RV ED volume / height	99R42			
RV_EDV_AI	RV ED volume / BSA	99R42			
RV_ESV	RV ES volume	99R42	8824-5	RV ES vol Imaging	LN
RV_ESV_ZValue	RV ES volume Z-Value	99R42			
RV_ESVI	RV ES volume / Height	99R42			
RV_ESV_AI	RV ES volume / BSA	99R42			
RV_SV	RV SV	99R42	8782-5	RV SV MRI	LN
RV_SV_ZValue	RV SV Z-Value	99R42			
RV_SVI	RV SV / height	99R42			
RV_SV_AI	RV SV / BSA	99R42	8802-1	RV SVI MRI	LN
RV_EF	RV EF	99R42	8818-7	RV EF MRI	LN
RV_EF_ZValue	RV EF Z-Value	99R42			
RV_CardiacOutput	RV Cardiac Output	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
RV_CardiacOutput_ZValue	RV Cardiac Output Z-Value	99R42			
RV_CardiacIndex	RV Cardiac Index	99R42	8762-7	RV Cardiac index MRI	LN
RV_Mass	RV Mass	99R42			
RV_Mass_ZValue	RV Mass Z-Value	99R42			
RV_MI	RV Mass / height	99R42			
RV_Mass_AI	RV Mass / BSA	99R42			
RV_MyoMass_Diast	RV Myo Mass Diastole	99R42			
RV_MyoMass_Diast_I	RV Myo Mass Diastole / height	99R42			
RV_MyoMass_Diast_AI	RV Myo Mass Diastole / BSA	99R42			
RV_MyoMass_Syst	RV Myo Mass Systole	99R42			
RV_MyoMass_Syst_I	RV Myo Mass Systole / height	99R42			
RV_MyoMass_Syst_AI	RV Myo Mass Systole / BSA	99R42			
RV_HeartRate	RV Heart Rate	99R42			
TAPSE	TAPSE	99R42			
ARVC/D (original)					
SevRVDilatationNormLV	Severe RV Dilatation and Reduced RV-EF Norm LV	99R42			
SevRVDilatationNormLV_LocInRV	SevRVDilatationNormLV Location in RV	99R42			
LocalizedRVAneurysms	Localized RV aneurysms with early diastolic bulge	99R42			
LocalizedRVAneurysms_LocInRV	LocalizedRVAneurysms Location in RV	99R42			
SevSegRVDilatation	Severe segmental RV dilatation	99R42			
SevSegRVDilatation_LocInRV	SevSegRVDilatation Location in RV	99R42			
MildGlobalRVDilatation	Mild global RV dilatation and/or reduced RV-EF norm LV	99R42			
MildGlobalRVDilatation_LocInRV	MildGlobalRVDilatation Location in RV	99R42			
MildSegRVDilatation	Mild segmental RV dilatation	99R42			
MildSegRVDilatation_LocInRV	MildSegRVDilatation Location in RV	99R42			
RegionalRVHypokinesis	Regional RV hypokinesis	99R42			
RegionalRVHypokinesis_LocInRV	RegionalRVHypokinesis Location in RV	99R42			
RVTrabecularHypertrophy	RV trabecular hypertrophy	99R42			
FibroFattyReplacement	Fibro-fatty replacement	99R42			
ARVC/D (revised)					
REV2_RegRVAkinesia	Regional RV akinesia	99R42			
REV2_RegRVDyskinesia	Regional RV dyskinesia	99R42			
REV2_DyssynchRVContraction	Dyssynchronous RV contraction	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
REV2_ModSevRVDilatation	Moderate Or Severe RV Dilatation	99R42			
REV2_MildRVDilatation	Mild RV Dilatation	99R42			
REV2_RVEjecFracMajor	RV ejection fraction <= 40%	99R42			
REV2_RVEjecFracMinor	40% < RV ejection fraction <= 45%	99R42			
REV2_MajorCriteria	ARVC/D Major Criteria	99R42			
REV2_MinorCriteria	ARVC/D Minor Criteria	99R42			
REV2_Microaneurysm	Microaneurysm	99R42			
REV2_SevSegRVDilatation	Severe segmental RV dilatation	99R42			
REV2_MildSegRVDilatation	Mild segmental RV dilatation	99R42			
REV2_AccordionSign	Accordion Sign	99R42			
REV2_FibfatReplacement	Fibro-fatty replacement	99R42			
REV2_TrabHypertrophy	RV trabecular hypertrophy	99R42			
Tissue Characterisation					
T2Ratio	T2 Signal Intensity ratio	99R42			
EarlyEnhancementRatio	Early Gadolinium Enhancement Ratio	99R42			
LateEnhancementMass	Late Gadolinium Enhancement Mass	99R42			
LateEnhLVMassRatio	Late Enhancement / LV Mass	99R42			
GreyzoneVolume	Greyzone Volume	99R42			
SalvagedAreaAtRisk	Salvaged Area At Risk	99R42			
MyoVolume	Myocardial Volume	99R42			
MyoMass	Myocardial Mass	99R42			
TotalEnhancedVolume	Total Enhanced Volume	99R42			
NoReflowVolume	No Reflow Volume	99R42			
NoReflowMass	No Reflow Mass	99R42			
ManuallyExcludedVolume	Manually Excluded Volume	99R42			
ManuallyExcludedMass	Manually Excluded Mass	99R42			
GreyzoneMass	Greyzone Mass	99R42			
ManuallyExcludedGreyzoneVolume	Manually Excluded Greyzone Volume	99R42			
ManuallyExcludedGreyzoneMass	Manually Excluded Mass	99R42			
GreyzoneRel	Greyzone Volume (rel)	99R42			
Valves and Flow (2D)					
ValveName	Valve Name	99R42			
SeriesID	Series ID	99R42			
ValveFunction	Valve Function	99R42			
BackgroundCorrectionMethod	Background Correction	99R42			
EncodingType	Encoding Type	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
ValveFlowHeartRate	ValveFlowHeartRate	99R42			
ValveTotalForwardFlow	Valve Total Forward Flow	99R42			
ValveForwardFlow	Valve Forward Flow	99R42			
ValveBackwardFlow	Valve Backward Flow	99R42			
ValvePeakVelocity	Valve Peak Velocity	99R42			
ValveRegurgitantFraction	Valve Regurgitant Fraction	99R42			
ValveOrificeArea	Valve Orifice Area	99R42			
ValveGradient	Valve Gradient	99R42			
MeanPressureGradient	Mean Pressure Gradient	99R42			
RawVelocityExtremaMax	Maximum Velocity	99R42			
RawVelocityExtremaMin	Minimum Velocity	99R42			
MinVol	Forward Flow	99R42			
MinVolEffective	Effective Forward Flow	99R42			
MaxAcceleration	Maximum Acceleration	99R42			
MinAcceleration	Minimum Acceleration	99R42			
RawFlowRateExtremaMax	Maximum Flow	99R42			
RawFlowRateExtremaMin	Minimum Flow	99R42			
NetPositiveVolume	Net Positive Volume	99R42			
NetNegativeVolume	Net Negative Volume	99R42			
COPulmonaryValve	Cardiac Output Pulmonary Valve	99R42			
COAorticValve	Cardiac Output Aortic Valve	99R42			
QPQS_Ratio	QP/QS Ratio	99R42	8743-7	Qp/Qs Imaging	LN
QPQSVolume	Qp/Qs Volume	99R42			
ShuntVolume	Shunt Volume	99R42			
ShuntVolumeMin	Shunt Vol/min	99R42			
QpSeriesID	Qp Series ID	99R42			
QsSeriesID	Qs Series ID	99R42			
Aorta					
AorticDiameter_Root	Aortic Root Diameter	99R42			
AorticDiameter_Ascending	Aortic Ascending Diameter	99R42			
AorticDiameter_Arch	Aortic Arch Diameter	99R42			
AorticDiameter_Descending	Aortic Descending Diameter	99R42			
AorticDiameter_AnnulusArea	Aortic Annulus Area	99R42			
AorticDiameter_SinusArea	Aortic Sinus Area	99R42			
Calcium Scoring					
LM_Volume	LM Ca Score Volume	99R42			
LM_Mass	LM Ca Score Mass	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
LM_Score	LM Ca Score	99R42			
LAD_Volume	LA Ca Score Volume	99R42			
LAD_Mass	LA Ca Score Mass	99R42			
LAD_Score	LA Ca Score	99R42			
CX_Volume	CX Ca Score Volume	99R42			
CX_Mass	CX Ca Score Mass	99R42			
CX_Score	CX Ca Score	99R42			
RCA_Volume	RC Ca Score Volume	99R42			
RCA_Mass	RC Ca Score Mass	99R42			
RCA_Score	RC Ca Score	99R42			
Total_Volume	Total Ca Score Volume	99R42			
Total_Mass	Total Ca Score Mass	99R42			
Total_Score	Total Ca Score	99R42			
CaThreshold	Ca Threshold	99R42			
CaMassCalibrationFactor	Ca Mass Calibration Factor	99R42			
SliceThickness	Ca Score Slice Thickness	99R42			
SliceDistance	Ca Score Slice Distance	99R42			
AgatstonScoreClassification	Agatston Score Classification	99R42			
CaScoreClassificationResult	Ca Score Classification Result	99R42			
T2*					
T2Star_ROI_SliceNum	T2Star ROI Slice Number	99R42			
T2Star_ROI_Name	T2Star ROI Name	99R42			
T2Star_ROI_T2Value	T2Star ROI T2 Value	99R42			
T2Star_ROI_T2Error	T2Star ROI T2 +/- Error	99R42			
T2Star_ROI_RSquare	T2Star ROI RSquare Value	99R42			
T2Star_ROI_FeValue	T2Star ROI Fe Value	99R42			
T2Star_ROI_FeError	T2Star ROI Fe +/- Error	99R42			
T2					
T2_ROI_SliceNum	T2 ROI Slice Number	99R42			
T2_ROI_Name	T2 ROI Name	99R42			
T2_ROI_T2Value	T2 ROI T2 Value	99R42			
T2_ROI_T2Error	T2 ROI T2 +/- Error	99R42			
T2_ROI_RSquare	T2 ROI RSquare Value	99R42			
T2_Global_Myo_Name	T2 Global Myo Name	99R42			
T2_Global_Myo_T2Value	T2 Global Myo T2 Value	99R42			
T2_Global_Myo_T2Error	T2 Global Myo T2 +/- SD	99R42			
T2_ROI_DCMManufacturer	T2 ROI Map Source	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
T2_ROI_T2MapSelected	T2 ROI Method	99R42			
T1					
T1_ROI_SliceNum	T1 ROI Slice Number	99R42			
T1_ROI_Series	T1 ROI Series	99R42			
T1_ROI_Name	T1 ROI Name	99R42			
T1_ROI_T1Value	T1 ROI T1 Value	99R42			
T1_ROI_T1Error	T1 ROI T1 +/- Error	99R42			
T1_ROI_RSquare	T1 ROI RSquare Value	99R42			
T1_Global_Myo_Series	T1 Global Myo Series	99R42			
T1_Global_Myo_Name	T1_Global_Myo_Name	99R42			
T1_Global_Myo_T1Value	T1_Global_Myo_T1 Value	99R42			
T1_Global_Myo_T1Error	T1 Global Myo T1 +/- SD	99R42			
T1_ROI_DCMManufacturer	T1 ROI Map Source	99R42			
T1_ROI_T1MapSelected	T1 ROI Method	99R42			
T1_ECV_Warning	ECV Warning	99R42			
Custom Measurements					
CustomMeasurementType	Custom Measurement Type	99R42			
CustomMeasurementLabel	Custom Measurement Label	99R42			
CustomMeasurementName	Custom Measurement Name	99R42			
CustomMeasurementValue	Custom Measurement Value	99R42			
CustomMeasurementError	Custom Measurement Error	99R42			
CustomMeasurementZValue	Custom Measurement Z-Value	99R42			
Findings and Summary					
Findings	Report Findings	99R42	11421-5	Phys CVS find	LN
NonNarrativeFindings	Non-Narrative Report Findings	99R42	11421-5	Phys CVS find	LN
TemplateFindings	Template Report Findings	99R42	11421-5	Phys CVS find	LN
Summary	Report Summary	99R42			
Cad_Rads	CAD-RADS	99R42			
Atrial Assessment					
Min_LA_Vol	Minimum Left Atrial volume	99R42			
Min_LA_Vol_I	Minimum Left Atrial volume / height	99R42			
Min_LA_Vol_AI	Minimum Left Atrial volume / BSA	99R42			
Max_LA_Vol	Maximum Left Atrial volume	99R42			
Max_LA_Vol_I	Maximum Left Atrial volume / height	99R42			
Max_LA_Vol_AI	Maximum Left Atrial volume / BSA	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
Min_RA_Vol	Minimum Right Atrial volume	99R42			
Min_RA_Vol_I	Minimum Right Atrial volume / height	99R42			
Min_RA_Vol_AI	Minimum Right Atrial volume / BSA	99R42			
Max_RA_Vol	Maximum Right Atrial volume	99R42			
Max_RA_Vol_I	Maximum Right Atrial volume / height	99R42			
Max_RA_Vol_AI	Maximum Right Atrial volume / BSA	99R42			
Min_LA_Area	Minimum Left Atrial area	99R42			
Max_LA_Area	Maximum Left Atrial area	99R42			
Min_RA_Area	Minimum Right Atrial area	99R42			
Max_RA_Area	Maximum Right Atrial area	99R42			
LeftAtrial_EF	Left Atrial Ejection Fraction	99R42			
Atrial_Method	Atrial Method	99R42			
RALVEDV	RA Vol at LVED volume	99R42			
RALVEDV_I	RA Vol at LVED volume / height	99R42			
RALVEDV_AI	RA Vol at LVED volume / BSA	99R42			
RALVESV	RA Vol at LVES volume	99R42			
RALVESV_I	RA Vol at LVES volume / height	99R42			
RALVESV_AI	RA Vol at LVES volume / BSA	99R42			
RARVEDV	RA Vol at RVED volume	99R42			
RARVEDV_I	RA Vol at RVED volume / height	99R42			
RARVEDV_AI	RA Vol at RVED volume / BSA	99R42			
RARVESV	RA Vol at RVES volume	99R42			
RARVESV_I	RA Vol at RVES volume / height	99R42			
RARVESV_AI	RA Vol at RVES volume / BSA	99R42			
RALVEDV_Area	RA Vol at LVED area	99R42			
RALVESV_Area	RA Vol at LVES area	99R42			
RARVEDV_Area	RA Vol at RVED area	99R42			
RARVESV_Area	RA Vol at RVES area	99R42			
LA_LVEDAV	LA Vol at LVED volume	99R42			
LA_LVEDAV_I	LA Vol at LVED volume / height	99R42			
LA_LVEDAV_AI	LA Vol at LVED volume / BSA	99R42			
LA_LVESAV	LA Vol at LVES volume	99R42			
LA_LVESAV_I	LA Vol at LVES volume / height	99R42			
LA_LVESAV_AI	LA Vol at LVES volume / BSA	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
LA_LVEDA_Area	LA Vol at LVED area	99R42			
LA_LVESA_Area	LA Vol at LVES area	99R42			
RA_LVEDA_Area	Right Atrial LV ED area	99R42			
RA_LVESA_Area	Right Atrial LV ES area	99R42			
RightAtrial_EF	Right Atrial Ejection Fraction	99R42			
LV_Atrial_Method	LV Atrial Method	99R42			
RA_Atrial_Method	RA Atrial Method	99R42			
RA_LV_Atrial_Method	RA LV Atrial Method	99R42			
LV_EDAV_Diameter	LA Diameter at LVED	99R42			
LV_ESAV_Diameter	LA Diameter at LVES	99R42			
RV_EDAV_Diameter	RA Diameter at LVED	99R42			
RV_ESAV_Diameter	RA Diameter at LVES	99R42			
Myocardial Strain					
SAX_Global_Circumferential_Strain_LV	SAX Peak Global Circumferential Strain - LV	99R42			
SAX_Global_Circumferential_Strain_RV	SAX Peak Global Circumferential Strain - RV	99R42			
LAX_Global_Longitudinal_Strain_LV	LAX Peak Global Longitudinal Strain - LV	99R42			
LAX_Global_Longitudinal_Strain_RV	LAX Peak Global Longitudinal Strain - RV	99R42			
Quantitative Perfusion					
RestMBF_LADTerritoryAverage	Rest MBF Territory Average - LAD	99R42			
RestMBF_RCATerritoryAverage	Rest MBF Territory Average - RCA	99R42			
RestMBF_LCxTerritoryAverage	Rest MBF Territory Average - LCx	99R42			
StressMBF_LADTerritoryAverage	Stress MBF Territory Average - LAD	99R42			
StressMBF_RCATerritoryAverage	Stress MBF Territory Average - RCA	99R42			
StressMBF_LCxTerritoryAverage	Stress MBF Territory Average - LCx	99R42			
RestRMBF_LADTerritoryAverage	Rest RMBF Territory Average - LA	99R42			
RestRMBF_RCATerritoryAverage	Rest RMBF Territory Average - RCA	99R42			
RestRMBF_LCxTerritoryAverage	Rest RMBF Territory Average - LCx	99R42			
StressRMBF_LADTerritoryAverage	Stress RMBF Territory Average - LAD	99R42			
StressRMBF_RCATerritoryAverage	Stress RMBF Territory Average - RCA	99R42			
StressRMBF_LCxTerritoryAverage	Stress RMBF Territory Average - LCx	99R42			
MPR_LADTerritoryAverage	MPR Territory Average - LAD	99R42			
MPR_RCATerritoryAverage	MPR Territory Average - RCA	99R42			
MPR_LCxTerritoryAverage	MPR Territory Average - LCx	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
rMPR_LADTerritoryAverage	"rMPR Territory Average - LAD	99R42			
rMPR_RCATerritoryAverage	"rMPR Territory Average - RCA	99R42			
rMPR_LCxTerritoryAverage	"rMPR Territory Average - LCx	99R42			
QP_StressHeartRate	Stress Heart Rate (QP)	99R42			
QP_RestHeartRate	Rest Heart Rate (QP)	99R42			
Stenosis					
Stenosis_Vessel	Stenosis Vessel	99R42			
Stenosis_Diameter	Stenosis Diameter	99R42			
Stenosis_Area	Stenosis Area	99R42			
Proximal_Diameter	Proximal Diameter	99R42			
Proximal_Area	Proximal Area	99R42			
Max_Stenosis_Diameter	Stenosis Diameter	99R42			
Max_Stenosis_Area	Stenosis Area	99R42			
Distal_Diameter	Distal Diameter	99R42			
Distal_Area	Distal Area	99R42			
Stenosis_Length	Stenosis Length	99R42			
Percent_Diameter	Percent Diameter	99R42			
Percent_Area	Percent Area	99R42			
Total_Plaque_Volume	Total Plaque Volume	99R42			
Total_Plaque_Burden	Total Plaque Burden	99R42			
Cp_Volume	Cp Volume	99R42			
Cp_Burden	Cp Burden	99R42			
Ncp_Volume	Ncp Volume	99R42			
Ncp_Burden	Ncp Burden	99R42			
Low_Density_Ncp_Volume	Low Density Ncp Volume	99R42			
Low_Density_Ncp_Burden	Low Density Ncp Burden	99R42			
Remodeling_Index	Remodeling Index	99R42			
Plaque_Ldncp_Threshold	Plaque Ldncp Threshold	99R42			
Stenosis_Phase	Stenosis Phase	99R42			
Stenosis V3					
Vessel_Name	Vessel Name	99R42			
Vessel_Phase	Vessel Phase	99R42			
Stenosis_Minimal_RCA	Minimal RCA (1-24%)	99R42			
Stenosis_Minimal_LAD	Minimal LM+LAD (1-24%)	99R42			
Stenosis_Minimal_LCX	Minimal LCX (1-24%)	99R42			
Stenosis_Minimal_Total	Minimal Total (1-24%)	99R42			
Stenosis_Mild_RCA	Mild RCA (25-49%)	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
Stenosis_Mild_LAD	Mild LM+LAD (25-49%)	99R42			
Stenosis_Mild_LCX	Mild LCX (25-49%)	99R42			
Stenosis_Mild_Total	Mild Total (25-49%)	99R42			
Stenosis_Moderate_RCA	Moderate RCA (50-69%)	99R42			
Stenosis_Moderate_LAD	Moderate LM+LAD (50-69%)	99R42			
Stenosis_Moderate_LCX	Moderate LCX (50-69%)	99R42			
Stenosis_Moderate_Total	Moderate Total (50-69%)	99R42			
Stenosis_Severe_RCA	Severe RCA (70-99%)	99R42			
Stenosis_Severe_LAD	Severe LM+LAD (70-99%)	99R42			
Stenosis_Severe_LCX	Severe LCX (70-99%)	99R42			
Stenosis_Severe_Total	Severe Total (70-99%)	99R42			
Stenosis_Occluded_RCA	Occluded RCA	99R42			
Stenosis_Occluded_LAD	Occluded LM+LAD	99R42			
Stenosis_Occluded_LCX	Occluded LCX	99R42			
Stenosis_Occluded_Total	Occluded Total	99R42			
Stenosis_CpVolume_RCA	Calcified Plaque RCA	99R42			
Stenosis_CpVolume_LAD	Calcified Plaque LM+LAD	99R42			
Stenosis_CpVolume_LCX	Calcified Plaque LCX	99R42			
Stenosis_CpVolume_Total	Calcified Plaque Total	99R42			
Stenosis_NcpVolume_RCA	Non-Calcified Plaque RCA	99R42			
Stenosis_NcpVolume_LAD	Non-Calcified Plaque LM+LAD	99R42			
Stenosis_NcpVolume_LCX	Non-Calcified Plaque LCX	99R42			
Stenosis_NcpVolume_Total	Non-Calcified Plaque Total	99R42			
Stenosis_LowDensityNcpVolume_RCA	(Non-Calcified Plaque) Low-Attenuation Plaque RCA	99R42			
Stenosis_LowDensityNcpVolume_LAD	(Non-Calcified Plaque) Low-Attenuation Plaque LM+LAD	99R42			
Stenosis_LowDensityNcpVolume_LCX	(Non-Calcified Plaque) Low-Attenuation Plaque LCX	99R42			
Stenosis_LowDensityNcpVolume_Total	(Non-Calcified Plaque) Low-Attenuation Plaque Total	99R42			
Stenosis_TotalPlaqueVolume_RCA	Total Plaque Volume RCA	99R42			
Stenosis_TotalPlaqueVolume_LAD	Total Plaque Volume LM+LAD	99R42			
Stenosis_TotalPlaqueVolume_LCX	Total Plaque Volume LCX	99R42			
Stenosis_TotalPlaqueVolume_Total	Total Plaque Volume Total	99R42			
Stenosis_TotalPlaqueBurden_RCA	Total Plaque Burden RCA	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
Stenosis_TotalPlaqueBurden_LAD	Total Plaque Burden LM+LAD	99R42			
Stenosis_TotalPlaqueBurden_LCX	Total Plaque Burden LCX	99R42			
Stenosis_TotalPlaqueBurden_Total	Total Plaque Burden Total	99R42			
StenosisVessel_TotalPlaqueVolume	Total Plaque Volume	99R42			
StenosisVessel_CpVolume	Calcified Plaque Volume	99R42			
StenosisVessel_NcpVolume	Non-Calcified Plaque Volume	99R42			
StenosisVessel_LowDensityNcpVolume	Low-Attenuation Plaque Volume	99R42			
StenosisVessel_TotalPlaqueBurden	Total Plaque Burden	99R42			
StenosisVessel_CpBurden	Calcified Plaque Burden	99R42			
StenosisVessel_NcpBurden	Non-Calcified Plaque Burden	99R42			
StenosisVessel_LowDensityNcpBurden	Low-Attenuation Plaque Burden	99R42			
StenosisVessel_MaxStenosis_Area_Name	Area Stenosis (Most Significant) Name	99R42			
StenosisVessel_MaxStenosis_Area	Area Stenosis (Most Significant)	99R42			
StenosisVessel_MaxStenosis_Diameter_Name	Diameter Stenosis (Most Significant) Name	99R42			
StenosisVessel_MaxStenosis_Diameter	Diameter Stenosis (Most Significant)	99R42			
StenosisVessel_HighestRemodelingIndex_Name	Remodeling Index (Highest) Name	99R42			
StenosisVessel_HighestRemodelingIndex	Remodeling Index (Highest)	99R42			
StenosisVessel_Name	Stenosis Vessel Name	99R42			
StenosisVessel_PercentArea	Percent Area	99R42			
StenosisVessel_PercentDiameter	Percent Diameter	99R42			
StenosisVessel_PlaqueThreshold_LAP	LAP	99R42			
StenosisVessel_PlaqueThreshold_NCP	NCP	99R42			
StenosisVessel_PlaqueThreshold_CP	CP	99R42			
StenosisVessel_PlaqueWarning	Plaque Warning	99R42			
Vascular CT					
VascularCT_MeasurementProtocol	Measurement Protocol	99R42			
VascularCT_MeasurementCategory	Measurement Category	99R42			
VascularCT_MeasurementName	Measurement Name	99R42			
VascularCT_MeasurementValue	Measurement Value	99R42			
VascularCT_MeasurementMax	Measurement Max	99R42			
VascularCT_MeasurementMin	Measurement Min	99R42			
VascularCT_MeasurementAvg	Measurement Avg	99R42			
VascularCT_MeasurementCRA	Measurement CRA	99R42			
VascularCT_MeasurementRAO	Measurement RAO	99R42			
VascularCT_VirtualGraftManufacturer	Virtual Graft Manufacturer	99R42			

Identifier OBX-3.1	Text OBX-3.2	Name of Coding System OBX-3.3	Alternate Identifier OBX-3.4	Alternate Text OBX-3.5	Name of Alternate Coding System OBX-3.6
VascularCT_VirtualGraftParts	Virtual Graft Parts	99R42			
VascularCT_VirtualGraftGraftTwist	Virtual Graft Graft Twist	99R42			
VascularCT_VirtualGraftAorticStartSlice	Virtual Graft Aortic Start Slice	99R42			
VascularCT_VirtualGraftInsertionSide	Virtual Graft Insertion Side	99R42			
VascularCT_VirtualDimensions	Virtual Graft Dimensions	99R42			
VascularCT_Comments	Vascular CT Comments	99R42			
Valves and Flow (4D)					
VF4D_Segmentation	Segmentation	99R42			
VF4D_RoiName	ROI Name	99R42			
VF4D_SeriesID	Series ID	99R42			
VF4D_OffsetCorrection	Offset Correction	99R42			
VF4D_AntiAliasing	Anti-Aliasing	99R42			
VF4D_ForwardFlow	Forward Volume	99R42			
VF4D_BackwardFlow	Backward Volume	99R42			
VF4D_TotalForwardFlow	Total Forward	99R42			
VF4D_NetPositiveVolume	Net Positive Volume	99R42			
VF4D_NetNegativeVolume	Net Negative Volume	99R42			
VF4D_RegurgitantFraction	Regurgitation Fraction	99R42			
VF4D_PeakVelocity	Peak Velocity	99R42			
VF4D_MinVol	Vol/min	99R42			
VF4D_MinVolEffective	Vol/min (Effective)	99R42			
VF4D_RawFlowRateExtremaMax	Maximum Flow	99R42			
VF4D_Gradient	Max Pressure Gradient	99R42			
VF4D_OrificeArea	Orifice Area	99R42			
VF4D_Function	Valve Assessment	99R42			
VF4D_Qp	Qp	99R42			
VF4D_Qs	Qs	99R42			
VF4D_QPQSVolume	Qp/Qs Volume	99R42			