



Molift Quick Raiser 205

A sit to stand hoist with
outstanding lifting capacity





Molift Quick Raiser 205

A sit to stand hoist with outstanding lifting capacity

Stable and compact

Molift Quick Raiser 205 is a state of the art sit-to-stand hoist, with an electrically adjustable leg base and an outstanding lifting capacity. The standing hoist is comfortable, stable and safe with exceptional manoeuvrability.

Quick Raiser 205 is compact, stable and exceptionally easy to manoeuvre. The unique wheel position, angled push bar and the extraordinary weight distribution of the base provides a small turning radius and easy manoeuvrability even in very narrow spaces. With the small castor option (50 mm), it is only 7,5 cm high from the floor, and fits under most beds and chairs while the electrically adjustable legs enables accessibility almost everywhere. The fixed-leg option benefits a very small footprint.

Natural pattern of movement

The Active Lifting Arm (four point sling bar) and the RgoSling Active provides extra support when hoisting from sit to stand. The inclined column encourages a natural pattern of movement and fulfils the hoisting needs of users with some weight bearing capacity. The Active Lifting Arm combined with the sling provides a comfortable, stable and safe hoisting solution.



Molift QuickRaiser 205

Model item no.

Battery	Variant	Adjustable leg	Fixed leg (FL)
SLA battery	Active, small casters	29000	29150
	StandUp, small casters	29001	29151
	Active, large casters	29100	29200
	StandUp	29101	29201
LI battery	Active, small casters	29002	29152
	StandUp, small casters	29003	29153
	Active, large casters	29102	29202
	StandUp, large casters	29103	29203

For complete variant matrix, contact your local distributor.

Weight

Maximum user weight (SWL)	205 kg
Total weight (FL)	44 (42,5) kg

Dimensions

External width of wheelbase	620 mm
Internal width of wheelbase	200 mm
Height of footboard:	
Castor size Ø 75/50mm	65 mm
Castor size Ø 100/75mm	90 mm

Technical

Battery SLA 24V	2,9 Ah
Charge time	6 hours
Rated performance	100 hoists/75 kg/500 mm
Lifting Range	600 mm
Turning radius	1150 mm
Hoisting speed 75 kg	35 mm/s
Hoisting speed 205 kg	30 mm/s



Clean the product with a solvent-free cleaning agent with a pH level between 5–9, or with a 70% disinfectant solution.

Features



Ergonomic design

The slanting column lifts the user up and forwards to imitate a natural movement pattern.



Comfort, safety and stability

The Active Lifting Arm (four point sling bar) and the RgoSling Active provides extra support when hoisting from sit to stand.



Automatic warning

Automatic warning for (annual) periodic inspection and service.



Easy to manoeuvre

Easy to operate in moving and handling situations; toilet, sit-to-stand and stand-up training.



The image includes Active lifting arm and RgoSling Active.


Accessories



StandUp lifting arm

Ideal to use with RgoSling Standup slings.

StandUp Lifting Arm 2910061

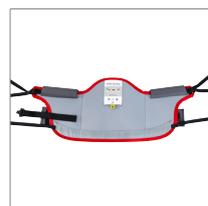

RgoSling Active & RgoSling StandUp

The slings are available in size XS to XL. For item numbers and measurements see www.etac.com



Footplate

This accessory benefits a 10 cm (4inch) elevated footrest. The attachment key is easy accessible and makes it convenient to lock or unlock the accessory to the Quick Raiser.

Footplate 2910065



Complete kit Bracket & Cushions

The alternative leg supports benefits an option of small or large cushions, and a leg strap for safety measurements

Kit with bracket and small cushions 2910059
Kit with bracket and large cushions 2910056



Set of inlay for footplate

Inlays with slip-resistant material to increase friction.

Set of inlay for footplate 2920151



Set of Cushions

Only contains cushions, no leg straps.

Set of small cushions 2910076
Set of large cushions 2910077



“Everyone should be able to live a free and independent life
and pursue their dreams regardless of any physical circumstances”

For the latest news and continuously updated
product information – visit: www.dlcana.com

DLCANA Inc.
323 Coventry Rd
Ottawa, ON K1K 3X6
Tel: 1-866-581-9030
Fax: 1-613-667-5655
customerservice@dlcana.com
www.dlcana.com

 **etac**[®]
Creating Possibilities