



# Favitor<sup>®</sup>

## DISPOSABLE RETRACTOR AND HOOKS

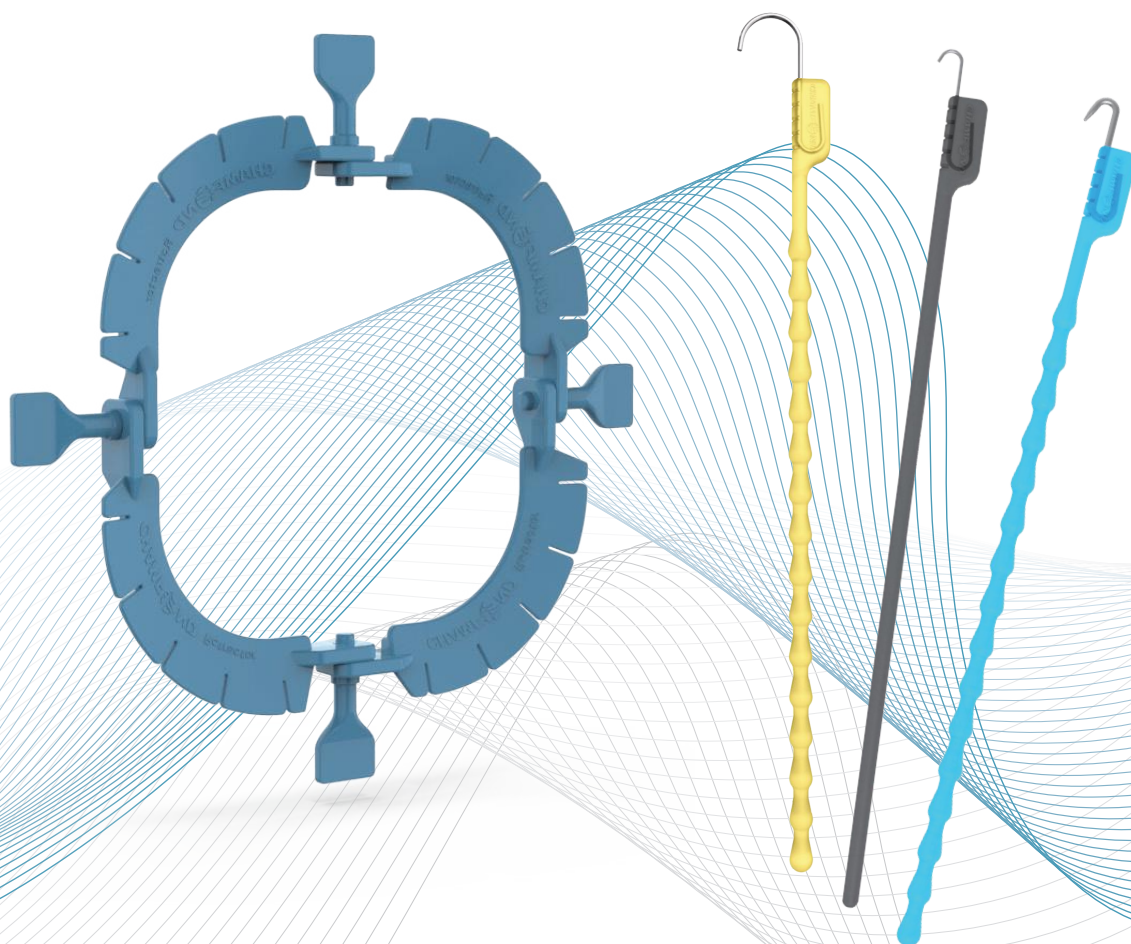


## Product model

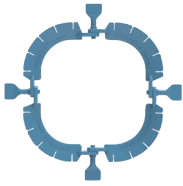


Model : C11  
Sizs : 14.4\*14.4cm  
Application : Colorectal surgery,  
general surgery, head and neck surgery




Model : C12  
Sizs : 16.6\*16.6cm  
Application : Gynecological surgery

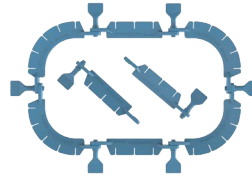

Model : C13  
Sizs : 25.4\*14.4cm  
Application : General surgery,  
neurosurgery, head and neck surgery

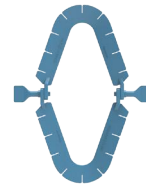

Model : C14  
Sizs : 18.6\*8.3cm  
Application : Head and neck surgery,  
otolaryngology, pediatric surgery

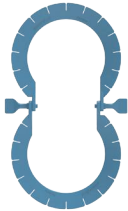

Model : C15  
Sizs : 25.3\*11.3cm  
Application : Gynecological surgery,  
head and neck surgery



Model : C16  
Sizs : 32.5\*18.3cm  
Application : Gynecological surgery,  
urological surgery



Model : C17  
Sizs : 22.0\*11.3cm  
Application : Gynecological surgery  
(pelvic stenosis/obesity)




Model : C18  
Sizs : 29.1\*18.3cm  
Application : Gynecological surgery,  
male urological surgery

## Product specification

Type	Description	Packaging	SKU No.
Single Sharp Hook	5mm Straight Hook	1 pcs/pack,50 packs/box	C21-1
		4 pcs/pack	C21-4
		8 pcs/pack	C21-8
	5mm Rain Drop Hook	1 pcs/pack,50 packs/box	C31-1
		4 pcs/pack	C31-4
		8 pcs/pack	C31-8
Single Blunt Hook	5mm Straight Hook	1 pcs/pack,50 packs/box	C22-05-1
		4 pcs/pack	C22-05-4
		8 pcs/pack	C22-05-8
	5mm Rain Drop Hook	1 pcs/pack,50 packs/box	C32-05-1
		4 pcs/pack	C32-05-4
		8 pcs/pack	C32-05-8
	12mm Straight Hook	1 pcs/pack,50 packs/box	C22-1
		4 pcs/pack	C22-4
		8 pcs/pack	C22-8
	12mm Rain Drop Hook	1 pcs/pack,50 packs/box	C32-1
		4 pcs/pack	C32-4
		8 pcs/pack	C32-8
Double Sharp Hook	7mm Straight Hook	2 pcs/pack	C23-2
		4 pcs/pack	C23-4
	7mm Rain Drop Hook	2 pcs/pack	C33-2
		4 pcs/pack	C33-4
Double Blunt Hook	12mm Straight Hook	2 pcs/pack	C24-2
		4 pcs/pack	C24-4
	12mm Rain Drop Hook	2 pcs/pack	C34-2
		4 pcs/pack	C34-4
Single Claw Rake Hook	16 mm x 20 mm Straight Hook	2 pcs/pack	C25-2
		4 pcs/pack	C25-4
	16 mm x 20 mm Rain Drop Hook	2 pcs/pack	C35-2
		4 pcs/pack	C35-4

# CLINICAL APPLICATION

## ● Anorectal Surgery

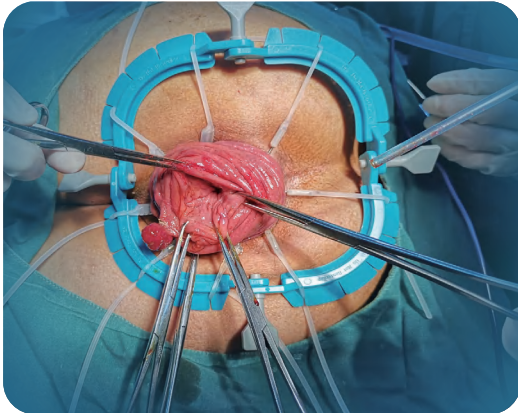


Recommended Model: C11+C31/C21




## Application


Middle and low rectal cancer  
Rectal mass  
Intestinal stromal tumor  
Rectal prolapse  
Anal fistula  
Hemorrhoids

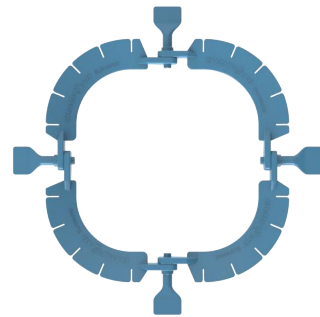

## ● ENT

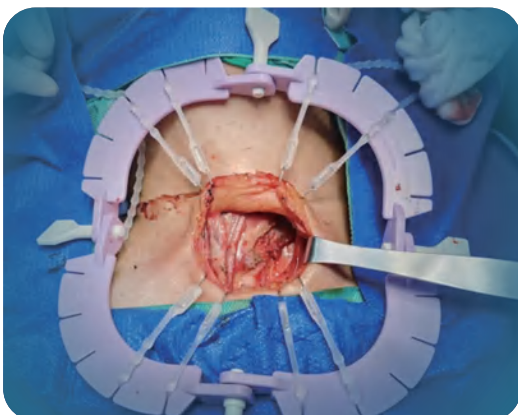


Recommended Model: C11/C13/C15+C21/C22/C25

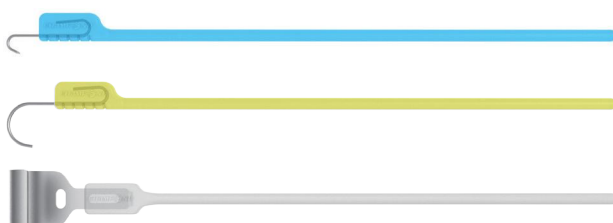


## Application

Breastectomy  
Parathyroid Surgery  
Tracheostomy / Ostomy  
Neck lymph node dissection  
Parotectomy  
Submandibular Glandectomy  
Lymph node biopsy

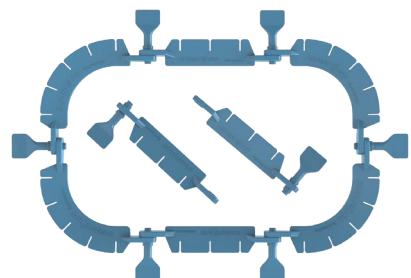

## ● Gynecology

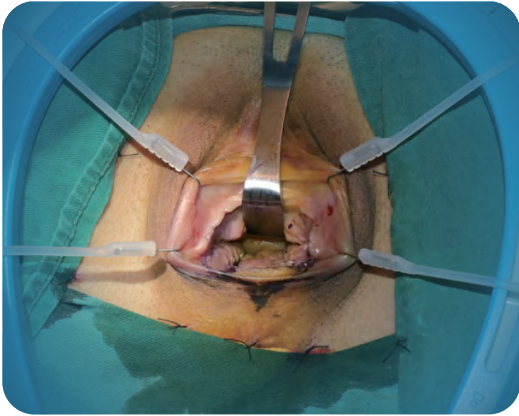


Recommended Model: C12/C16+C21



## Application

Rectovaginal Fistula  
Vaginal repair/reconstruction  
Pelvic Floor Repair/Reconstruction  
Vagina tightening

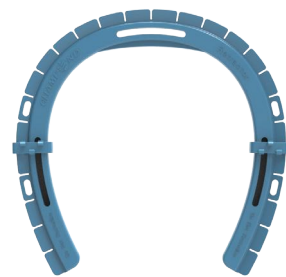

## ● Urology




Recommended Model: C16+C31/C32/C21/C22




## Application

Urethral stricture  
Urethrostomy

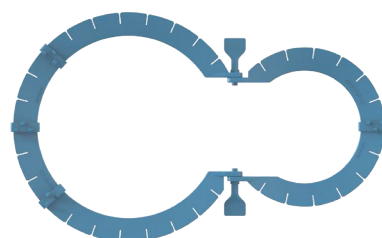




## ● Neurosurgery

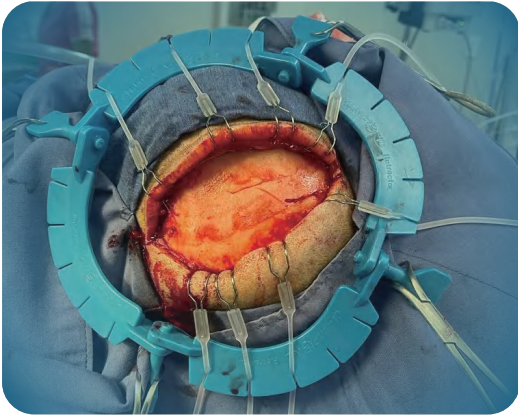


Recommended Model: C11+C31/C21/C22/C25

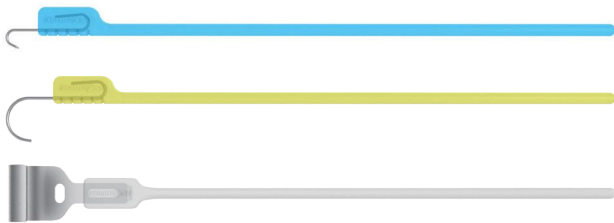


## Application

Convex meningioma  
Carotid endarterectomy  
Superficial Temporal Arteries Surgery  
Aneurysms  
Vestibulectomy

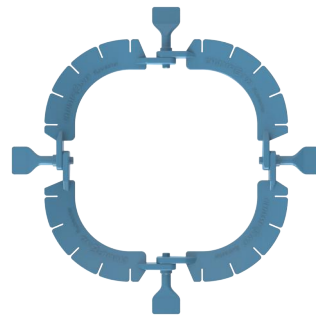

## ● Veterinary

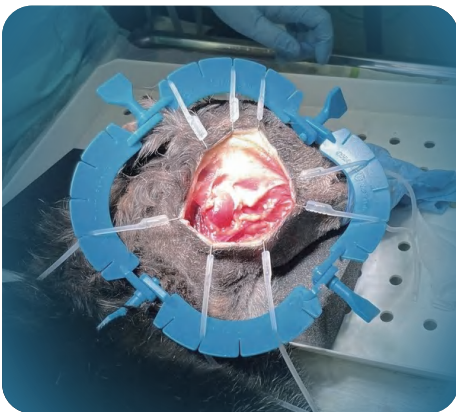

Recommended Model: C11+C21/C31




## Application

Cataract removal(canine)  
Liguinal hernia repair(canine)  
Laminectomy(canine)  
Microsurgery(rodent)  
Ovarian cyst(canine)  
Rumenotomy(ovine)

# PATENTED TECHNOLOGY

- Limit slot for joint part

Authorized announcement No. CN 210631251 U

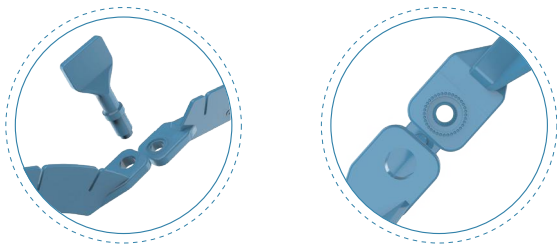


Advantage

Prevents deformation of retractor ring under great force  
Provides reliable and stable surgical site retraction

- Raindrop chain design


Authorized announcement No. CN 211155975 U

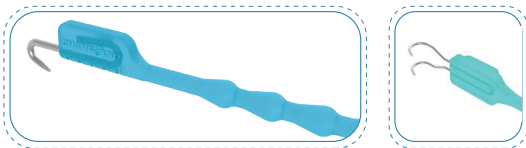

Advantage


Withstand greater pulling force, 20% more than the linear type  
Avoid sliding, maintain tissue retraction stability

- Catheter clips for urology

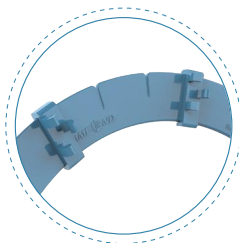


Easy to grip catheter, Improve the convenience of operation

- Velcro

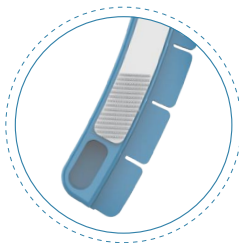

Attached on sterile surgical drapes, Improve stability on site

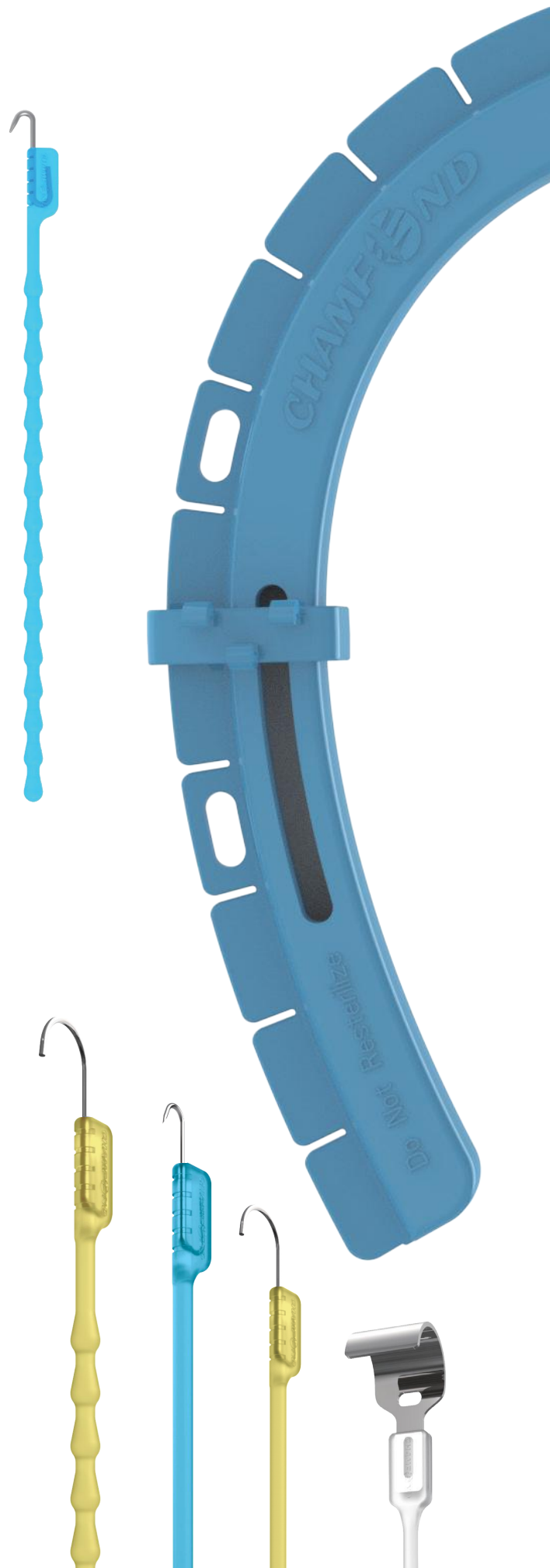

---

## DISPOSABLE RETRACTOR AND HOOKS

- Easy to deploy, surgical field of view exposed layer by layer
- 360° exposure of the surgical field of view, reducing space occupation, and accurate visualization of the entire operation
- Full automatic pulling, reducing assistant labor
- Flexible combination of different hooks and retractor rings, suitable for various of surgical operations
- Multi-point traction to avoid tissue strain or necrosis
- Single use to avoid cross-infection

---

## **FLEXIBLE COMBINATION FULLY EXPOSED**

---

Better Vision Better Care

**Chamfond Biotech Co., Ltd**

Tel:400-997-6833 Fax:025 5886 8783

E-mail:info@chamfondbiotech.com

**[www.chamfondbiotech.com](http://www.chamfondbiotech.com)**

A series of thin, light blue wavy lines that sweep across the bottom right corner of the page, creating a sense of motion and modern design.