



Molift Mover 180

A lightweight and compact hoist, suitable in home care and institutional environments

Easy to manoeuvre

The Molift Mover 180 features a lightweight aluminium construction, making it easy to manoeuvre. Its large hoisting interval of 27–163 cm allows for hoisting from both high surfaces and the floor. Both 2-point and 4-point sling bars are available and can easily be changed without tools thanks to the quick release coupling,

Equipped with 75 mm castors, it provides excellent manoeuvrability in tight spaces. For movement over uneven surfaces, 100 mm castors are available, enhancing stability and ease of use.



molift[®]



Hand control

Robust and washable, Indicators for battery, overload, and service.



Unique push bar

The push bar ensures ergonomic, natural grip with multiple options.



Sling bar options

2-point and a 4-point sling bars are available in several sizes.



Easy to manoeuvre

Lightweight and compact design makes it easy to manoeuvre.

Features

Ergonomic handling

The Molift Mover 180 is equipped with Molift's unique push bar that provides multiple and flexible grip possibilities for improved caregiver ergonomics.

Life expectancy

The hoist has an expected lifetime of 30,000 cycles or 15 years, whichever comes first.

Battery

The battery features advanced Li-Ion technology and a safe, sealed design makes it washable (IPX6). It weighs only 0.9 kg, and can be replaced with a single hand.

The battery can be charged directly on the hoist, and an external charger is also available. The charging time is only 3.5 hours. Fully charged it delivers 160 hoists at 75 kg per, and is compatible with all Molift mobile hoists.

With a life expectancy of 1,700 charging cycles and a 5-year warranty, the battery is built to last.



Speed boost

The speed boost feature enables increased speed, when the actuator is unloaded.

Hand control


The hand control is robust and well-proven. It features a battery level indicator, overload indicator, and service indicator to ensure optimal performance and safety.

Designed for durability, the cable can be easily replaced if damaged, and the hand control is washable (IPX6).

Molift Service Tool

Molift Mover 180, like all Molift hoists, is equipped with a unique service notification system that notifies when service is required. The integrated service software can also display data about hoist usage and can assist service technicians in troubleshooting.

Additionally, the hoist is equipped with Bluetooth connectivity, enabling wireless connection to Molift Service Tool.

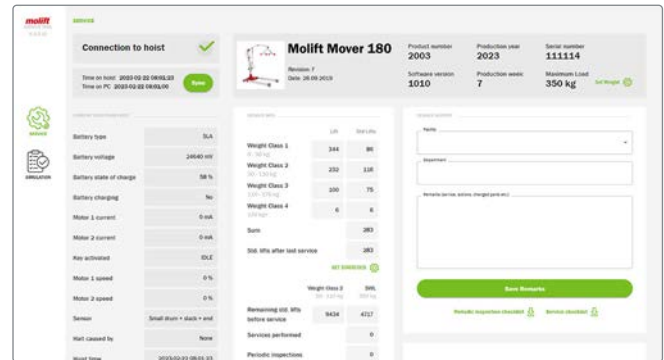

Product	Item no.	Dimensions	Technical
Molift Mover 180	M1801	Outside measurement with closed legs	670 mm
Molift Mover 180	M1802	Outside measurement with open legs	1145 mm
- excluding sling bar (ordered separately)		Inside measurement with closed legs	555 mm
		Inside measurement with open legs	1025 mm
Material		Technical	
Aluminium, steel, plastic		Battery	Li-Ion 25.2V, 2.85Ah
Weight		Battery charge time	3.5 hours
Safe working load (SWL)	180 kg	Hoisting speed	46 mm/sec
Weight of unit (excl. battery)	31 kg	Rated performance	160 lifts/75 kg
Dimensions		Cleaning	
Height	1940 mm		
Length	1260 mm	Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution	
Leg height	105 mm		
Castor size	75/100 mm		
Hoisting interval	1365 mm		
Hoisting min-max with 2-point sling bar	360 - 1730 mm		
Hoisting min-max with 4-point sling bar	265 - 1630 mm		

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com