




Tel: 1-613-258-3344  
Email: [inquires@tlcmedicalsupply.ca](mailto:inquires@tlcmedicalsupply.ca) & [mario@tlcmedicalsupply.ca](mailto:mario@tlcmedicalsupply.ca)

## LORI / LORI ELECTRIC

Lori stander enables to verticalize the patient both in the prone and supine position.

The device offers the anatomically-correct abduction of lower limbs with angle measurement both for the left and the right leg ( $0^{\circ}$ –  $35^{\circ}$ ). This function ensures correct placement of the hip joint and hence we can prevent hip dysplasia. The device allows daily weight bearing exercises.

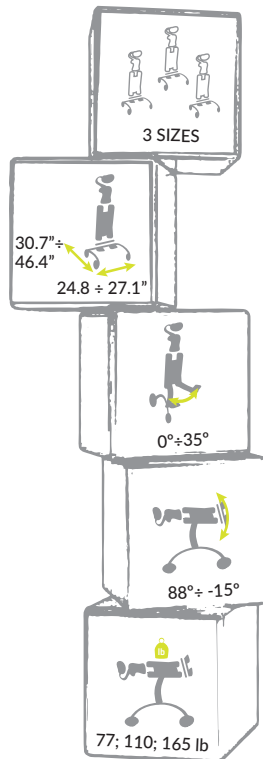

**COLOURS**



**IMPORTANT FEATURES**



A height-adjustable trunk and hip support system together with lateral supports secure and position the patient in safe and stable therapeutical position.



3D back rotational support, patented system, that supports patient's back, hips and pelvic.



Legs are protected by knee supports, legs abduction system and 3D foot rotational platforms. All equipped with independent adjustment for left and right lower limb.



The tray with a bowl. Use it for play and therapy.



Lori Electric equipped with an electric actuator for easy daily use.



Headrest made of flexible material firmly holds the head in supine or prone standing position.

**TECHNICAL DATA**

Base width/length

**SIZE 1. baby**

63/78 cm / 24.8"/30.7"

**SIZE 2. junior**

63/78 cm / 24.8"/30.7"

**Size 3. teen- Lori Electric**

69/118 cm / 27.1"/46.4"

Total height (without headrest)

94 cm / 37"

120 cm / 47.2"

160 cm / 63"

Range of verticalization

90°÷-8°

90°÷-15°

88°÷-2°

Maximum user weight

35 kg / 77 lb

50 kg / 110 lb

75 kg / 165 lb