

DEDICATED TO REVOLUTIONIZING THE POWER CHAIR EXPERIENCE



"THE EDGE 4 CARRIES YOU FORWARD WITH POWER,
GRACE, AND THE KIND OF CONFIDENCE THAT
MAKES THE WORLD FEEL WIDE OPEN."

— **STEWART LUNDY**
DENVER, CO



FDA Class II Medical Device



quantumrehab.com





FDA Class II
Medical Device



FDA Class II
Medical Device





W E L C O M E T O **QUANTUM**[®]

Quantum Rehab[®] was born out of the desire to delight customers with the most advanced, consumer-inspired complex rehab power wheelchairs and related technologies possible.

Clinically-focused and consumer-inspired, Quantum[®] products are engineered and manufactured to meet the widest range of complex rehab needs with functional, durable, stylized, high-performance designs to increase daily functionality and quality of life. At Quantum, consumer needs and wishes are the driving force.

We're dedicated to not just meeting medical and clinical needs, but also quality-of-life needs. From the most advanced power seating for pressure management to USB ports, Bluetooth[®] and fender lights, no consumer need is overlooked.

Headquartered in Duryea, Pennsylvania, Quantum also has operations in Australia, Canada, China, Italy, France, Germany, the Netherlands, New Zealand, and the United Kingdom.

quantumrehab.com



Quantum[®] Power Chairs are
FDA Class II Medical Devices
designed to aid individuals with
mobility impairments.

FDA Class II
Medical Device





Safer transfers



Improved reach



FDA Class II
Medical Device

"COOKING IS SO MUCH EASIER NOW THAT I CAN REACH THINGS. I CAN'T IMAGINE MY LIFE WITHOUT iLEVEL."

STEPHANIE WOODWARD
ROCHESTER, NY

iLEVEL®

POWER SEAT ELEVATION

Standing Height + Speed + Stability = Authentic iLevel (like no other)

iLevel® power seat elevation technology features a revolutionary new power chair design for a better experience. iLevel improves daily functionality and quality of life by allowing operation of the power chair with the seat elevated while at walking speed.

FEATURES

- 12" of power adjustable seat height while at walking speed; maximum of 10" height on the Q6 Edge® HD
- iLevel's patented Extra Stability Technology® enhances safety for transfers, reach, and numerous other activities of daily living
- Available on the Stretto, Edge® 4, Edge® 3, and Q6 Edge® 2.0
- Front and rear LED fender lights come standard with iLevel®



FDA Class II
Medical Device



FDA Class II
Medical Device



Social
interaction
at iLevel®

FDA Class II
Medical Device

FDA Class II
Medical Device

PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS



OPTIONAL
Quantum® backup camera

Q-Logic 3 controller with Bluetooth® and USB charger

OPTIONAL
TRU-Balance® 4 power positioning with anterior tilt, 12" iLevel® and memory seating

Stealth Products® cushion featuring exclusive Coolcore® technology

Switch-operated rear LED accent lights integrated into fenders

Power articulating foot platform

Independent Smooth Ride Suspension

Switch-operated LED accent lights integrated into the casters for increased light projection

New battery tray for improved rear battery access

7" front and rear casters

The Edge® 4 is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

New easy-access rear battery tray



Integrated LED lights



"THE EDGE 4 ISN'T JUST MOBILITY—IT'S MOMENTUM. IT CARRIES YOU FORWARD WITH POWER, GRACE, AND THE KIND OF CONFIDENCE THAT MAKES THE WORLD FEEL WIDE OPEN."

STEWART LUNDY
DENVER, CO

EDGE® 4

MID-WHEEL DRIVE

The Edge® 4 redefines full-size mid-wheel drive with 7" caster wheels, group 24 batteries, LED accent lights, and improved battery access. With standard 4-pole motors, the Edge 4 can reach a top speed of 7 mph, delivering exceptional speed and torque. SRS (Smooth Ride Suspension) offers enhanced comfort, and optional iLevel® provides up to 12" of seat elevation at a walking speed of 4.5 mph.

SPECS

MAXIMUM SPEED: 7 mph (4.5 mph when elevated with iLevel®)


WEIGHT CAPACITY: 300 lbs. (See Q6 Edge® HD for increased capacity)

TURNING RADIUS: 22.8"

POWER SEAT ELEVATION: Optional 12" iLevel

COLORS

* See website for a complete list of color options.

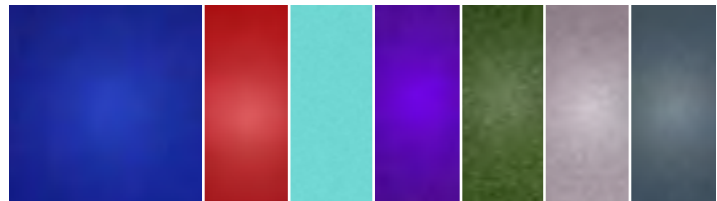

ROYALE BLUE

RED INFERNO

BERMUDA BLUE

COSMIC PURPLE

FOREST GREEN

ICE PINK

STEEL BLUE



Independent Smooth
Ride Suspension



PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS



OPTIONAL

TRU-Balance® 4 Power Positioning Systems provide a choice of 10, 20, or 30 degrees of power anterior tilt.

OPTIONAL

12" iLevel® power seat elevation at 4.5 mph (standard on TB4)

Featuring Extra Stability Technology® for safer transfers, reach and other activities of daily living



Q-Logic 3 controller with Bluetooth® and USB charger

Stealth Products® cushion featuring exclusive Coolcore® technology

Upgraded SRS (Smooth Ride Suspension)

The Edge® 3 is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

USB charger (standard)



Front and rear LED fender lights for seeing and being seen (standard)



"THE EDGE 3 WITH iLEVEL HAS RESTORED MY INDEPENDENCE AT BOTH WORK AND HOME. iLEVEL HAS ALLOWED ME TO GET MORE INVOLVED AND START A DIALOGUE WITH PEOPLE, SOMETHING I WOULDN'T HAVE DONE BEFORE."

MERLISHA HENDERSON
PHOENIX, AZ

EDGE® 3

MID-WHEEL DRIVE

The Edge® 3 offers upgraded SRS (Smooth Ride Suspension), which provides a smooth, comfortable ride and advanced stability. Optional iLevel® is available at 4.5 mph with up to 12" of power adjustable seat height.

SPECS

MAXIMUM SPEED: 6.25 mph (7 mph with optional 4-pole motors and Q-Logic 3 / 4.5 mph when elevated with iLevel®)

WEIGHT CAPACITY: 300 lbs. (See Q6 Edge® HD for increased capacity)

TURNING RADIUS: 20.5"

POWER SEAT ELEVATION: Optional 12" iLevel

COLORS

* See website for a complete list of color options.



GREEN MACHINE

RED CORVETTE

DEEP PURPLE

LEMON CRUSH

RASPBERRY BERET

EVOLUTION ORANGE

BLUE JEAN

AVAILABLE WITH iLEVEL
POWER SEAT ELEVATION

Drive wheel
color accents
(standard)



PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS

NARROWEST POWER
BASE IN AMERICA AT

20.47"

WIDTH WITH 12.5"
DRIVE WHEELS



OPTIONAL
TRU-Balance® 4 Power
Positioning Systems provide
a choice of 10, 20, or 30
degrees of power anterior tilt.

Q-Logic 3 controller
with Bluetooth® and
USB charger

OPTIONAL
12" iLevel® power seat
elevation at 3.5 mph
(standard on TB4)

Stealth Products® cushion
featuring exclusive Coolcore®
technology

Toolless removal of
rear shroud with
battery access door

Front and rear LED
fender lights

Independent SRS
(Smooth Ride
Suspension)



The Stretto is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

Independent Smooth
Ride Suspension

The Stretto is the ideal
pediatric power chair solution



"TIGHT ROOMS, NARROW DOORWAYS, OUT-OF-REACH CABINETS AND VAN RAMPS... NOW I CAN GO WHERE I WANT TO GO, ALL BECAUSE OF THE STRETTO!"

BRYSON FOSTER
CONCORD, NC

STRETTO

MID-WHEEL DRIVE

The Edge 3 Stretto® is the narrowest, most maneuverable power base in America! In addition to its narrow width, it is equipped with independent SRS (Smooth Ride Suspension), making it ideal for children, teenagers and small adults. The Stretto is available with 12.5" or 14" drive wheels and optional 12" of iLevel® power adjustable seat height.

SPECS

MAXIMUM SPEED: 6.25 mph (3.5 mph when elevated with iLevel®)

WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 21"

POWER SEAT ELEVATION: Optional 12" iLevel

COLORS

* See website for a complete list of color options.



BLUE JEAN

RED CORVETTE

DEEP PURPLE

LEMON CRUSH

RASPBERRY BERET

EVOLUTION ORANGE

GREEN MACHINE

AVAILABLE WITH iLEVEL
POWER SEAT ELEVATION

Toolless access
to electronics



PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS

THE EVOLUTION OF REAR WHEEL DRIVE



Q-Logic 3 Advanced Drive Controls with Bluetooth® and USB charger

OPTIONAL
Quantum backup camera

Stealth Products® cushion featuring exclusive Coolcore® technology

OPTIONAL
10" safe seat elevation at up to 3.5 mph

Independent SRS (Smooth Ride Suspension)

LED lighting in front caster arms and rear fenders

The R-TRAK is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

Independent SRS (Smooth Ride Suspension)

LED lighting in front caster arms and rear fenders



R-TRAK

As a global innovator, we are creating a product that offers confidence and stability elevated at the speed of life. Our state-of-the-art facility has a robust team of experienced technicians who custom build each R-TRAK specifically to the individual consumer's needs.



Scan to learn about the design and development of the R-TRAK.



R-TRAK

REAR-WHEEL DRIVE

The R-TRAK revolutionizes the rear-wheel drive power chair experience with independent SRS suspension on all four wheels, providing a smoother, more maneuverable ride. Standard features include fender lights, USB charger and drive wheel color accents. The R-TRAK is also available with our 10" safe seat elevation, allowing the client to drive up to 3.5 mph while elevated.

SPECS

MAXIMUM SPEED: Up to 7 mph (3.5 mph when elevated)


WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 25.7"

POWER SEAT ELEVATION: Optional 10" safe seat elevation

COLORS

* See website for a complete list of color options.

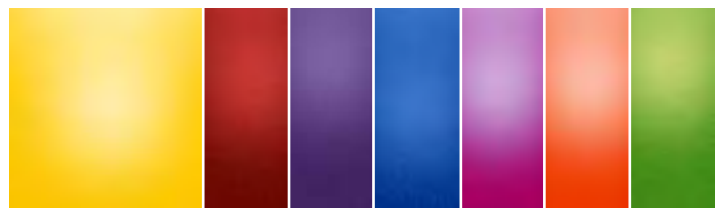

LEMON CRUSH

RED CORVETTE

DEEP PURPLE

BLUE JEAN

RASPBERRY BERET

EVOLUTION ORANGE

GREEN MACHINE

Optional Quantum backup camera can be mounted anywhere



PROUDLY MADE IN
AMERICA
OF FOREIGN AND DOMESTIC PARTS



450-LB. WEIGHT CAPACITY

available on the
4FRONT® 2 HD

Optional backup camera with 3.5" display can be mounted in multiple locations

OPTIONAL

TRU-Balance® 4 Power Positioning Systems provide a choice of 10°, 20°, or 30° of power anterior tilt.

OPTIONAL

10" power adjustable safe seat elevation at 3.5 mph (standard on TB4)

Redesigned caster arms for improved durability

STC (Smart Traction Control)



Q-Logic 3e controller featuring a magnetic charging port (Q-Logic 3 also available)

Stealth Products® cushion featuring exclusive Coolcore® technology

Front and rear LED fender lights for seeing and being seen (standard)

Automotive-grade SRS (Smooth Ride Suspension)

The 4Front® 2 is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

450-lb. weight capacity is available on the **4FRONT® 2 HD**

Automotive-grade SRS (Smooth Ride Suspension)



"RIDING IN THE 4FRONT 2 IS UNLIKE ANY CHAIR I HAVE EVER USED. THE 4FRONT 2 PUSHES THE LIMIT SET BY ITS PREDECESSOR BEYOND MY EXPECTATIONS."

CHRISTIAN BUDNEY
SCRANTON, PA

4FRONT® 2

FRONT-WHEEL DRIVE

Featuring STC (Smart Traction Control)

The 4Front® 2 is the next generation of the 4Front® Power Chair! New Smart Traction Control (STC), along with Smooth Ride Suspension (SRS), provide stability for smooth obstacle transitions and unmatched driving performance. Additional new features include 4-pole motors and redesigned caster arms for improved durability and aesthetics and reduced weight.

SPECS

MAXIMUM SPEED: 6.25 mph (3.5 mph when elevated)

WEIGHT CAPACITY: 300 lbs. (450 lbs. on 4Front® 2 HD)

TURNING RADIUS: 25.9"

POWER SEAT ELEVATION: Optional 10" safe seat elevation

COLORS

* See website for a complete list of color options.



EVOLUTION ORANGE

RED CORVETTE

DEEP PURPLE

LEMON CRUSH

RASPBERRY BERET

GREEN MACHINE

BLUE JEAN

AVAILABLE WITH 10"
SAFE SEAT ELEVATION

LED front and rear fender lights
(front: white / rear: red)





MAKE YOUR OWN PATH.



Available with TRU-Balance® 3 Power Positioning (shown) or Q-Captain Seating

Optional Quantum® backup camera

OPTIONAL 10" safe seat elevation at up to 3.2 mph

Rear LED lights



Independent SRS (Smooth Ride Suspension)

Q-Logic 3e controller with Bluetooth®

Stealth Products® cushion featuring Coolcore® technology

Front LED lights

15" pneumatic tires for superb control over rough terrain

The Quantum® Outback is an FDA Class II Medical Device designed to aid individuals with mobility impairments.

FRONT LED LIGHTING

REAR LED LIGHTING

Q-LOGIC 3e CONTROLLER WITH BLUETOOTH®



"THE OUTBACK LETS ME LIVE LIFE THE WAY I WANT. WHETHER I'M WORKING ON THE FARM OR GOING FOR A HIKE WITH MY FRIENDS, THE OUTBACK ALLOWS ME TO GO WHEREVER I WANT TO GO."

BRYAN ANDERSON
NASHVILLE, TN

QUANTUM® OUTBACK

4X4 ALL-TERRAIN

Quantum Rehab® welcomes adventurers and nature enthusiasts to unleash their wild side with the release of our cutting-edge 4x4 all-terrain power chair, the Quantum® Outback! Designed for those who dare to defy conventional limits, this rugged mobility solution promises to elevate any outdoor adventure, whether navigating wooded trails or exploring uncharted terrain.

SPECS


MAXIMUM SPEED: Up to 5 mph (6 mph with optional upgrade)

WEIGHT CAPACITY: 300 lbs. (350 lbs. with Q-Captain Seating)

TURNING RADIUS: 49.2"

POWER SEAT ELEVATION: Optional 10" safe seat elevation

COLORS

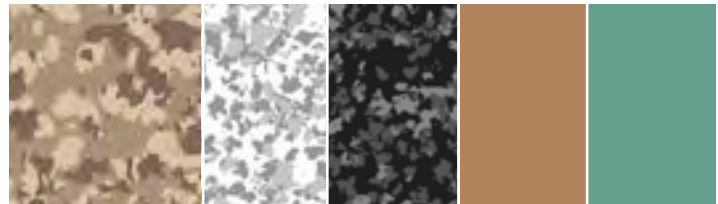

DESERT CAMO

WHITE CAMO


BLACK CAMO

TAN (MATTE)

CADET (MATTE)



INDEPENDENT SMOOTH RIDE SUSPENSION

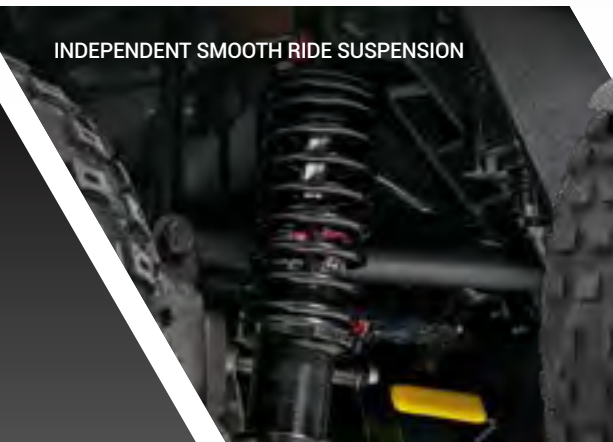

FDA Class II
Medical Device



Q6 EDGE® 2.0

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 6.25 mph (3.5 mph when elevated)


WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 20.5"

POWER SEAT ELEVATION: Optional 12" iLevel®

COLORS

* See website for a complete list.

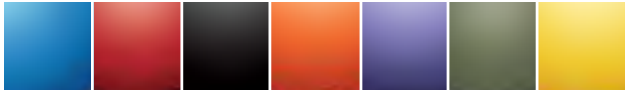

FDA Class II
Medical Device



Q6 EDGE® HD

BARIATRIC MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 4.7 mph (3.5 mph when elevated)

WEIGHT CAPACITY: 450 lbs. (400 lbs. with iLevel®)

TURNING RADIUS: 19.75"

POWER SEAT ELEVATION: Optional HD 10" iLevel®

COLORS

* See website for a complete list.



FDA Class II
Medical Device



FDA Class II
Medical Device

Q1450

BARIATRIC FRONT-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 4.5 mph

WEIGHT CAPACITY: 600 lbs.

TURNING RADIUS: 26.5"

COLORS



FDA Class II
Medical Device



J4E

MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 5 mph

WEIGHT CAPACITY: 300 lbs.

TURNING RADIUS: 20.5"

POWER SEAT ELEVATION: Optional 10" Safe Seat Elevation

COLORS



J4E
shown with
optional Q4
Power Tilt



FDA Class II
Medical Device

FDA Class II
Medical Device



J4E HD

BARIATRIC MID-WHEEL DRIVE

SPECS

MAXIMUM SPEED: 5 mph

WEIGHT CAPACITY: 400 lbs.

TURNING RADIUS: 20"

POWER SEAT ELEVATION: Optional 10" Safe Seat Elevation

COLORS



J4E HD
shown with
optional Q4
Power Tilt



FDA Class II
Medical Device

#1 FOR REHAB POWER

BUILT TOUGH

Quantum® bariatric power chairs all offer an innovative design with superior outdoor performance and superb indoor maneuverability. Our bariatric line of power chairs offers speeds as high as 5 mph and weight capacities up to 600 lbs.

All of our bariatric power chairs can be equipped with power positioning for greater comfort, such as TRU-Balance® HD Power Tilt. This seating system can assist with pressure relief. To provide additional support, Stealth® cushions and backs are available to provide lumbar support or coccyx and sacral relief. All Quantum bariatric power chairs are available with Q-Logic controls featuring Bluetooth® technology.



QUANTUM®

B A R I A T R I C S

Quantum® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.

Q6 Edge® HD

450 LBS.
WEIGHT CAPACITY



FDA Class II
Medical Device

Q1450

600 LBS.
WEIGHT CAPACITY



FDA Class II
Medical Device

| | Q6 Edge® HD | Q1450 | J4 HD | 4Front® 2 HD |
|------------------------|-------------|----------|-----------|--------------|
| Capacity | 450 lbs. | 600 lbs. | 400 lbs. | 450 lbs. |
| Top Speed | 4.7 mph | 4.5 mph | 5 mph | 5 mph |
| Turning Radius | 19.75" | 26.5" | 20" | 24.5" |
| Power Elevation | Yes (10") | N/A | Yes (10") | Yes (10") |

*** Please consult the back of the catalog for detailed specifications.*



J4E HD
400 LBS.
 WEIGHT CAPACITY

FDA Class II
 Medical Device



4Front® 2 HD
450 LBS.
 WEIGHT CAPACITY

FDA Class II
 Medical Device

TRU-Balance® 4 Power Positioning Systems

WITH ANTERIOR TILT AND MEMORY SEATING

TRU-Balance® 4 anterior tilt is a real game changer. Users will enjoy the ease in which they can perform ADLs such as accessing sinks and cooktop surfaces, meal preparation, and reaching items in cabinets.

FEATURES

- Maximum weight capacity of up to 300 lbs.
- Seat elevation combined with anterior tilt provides optimal access.
- Up to 8 customizable memory seating positions recall the saved positions of all actuators, offering the choice between sequential or synced movement for each seat function.
- Memory seating profiles can be labeled for easy recall of most used positions.



Chest Bar



FDA Class II
Medical Device

TRU-Balance® 4
Power Recline



TRU-Balance® 4
Power Posterior Tilt



DUAL ACTUATOR POWER ARTICULATING FOOT PLATFORM

Lowers to assist with transfers

Articulation range of 0-7.25"

Elevation range of 85-165 degrees

Only available on TRU-Balance® 4

TRU-BALANCE® 4

POWER POSITIONING SYSTEMS

FEATURES

- Anterior tilt options: 10, 20 and 30 degrees
- Up to eight customizable memory seating positions
- Chest bar and knee blocks available (knee blocks and chest support required with 20 or 30 degrees of anterior tilt)
- Available on the Stretto, 4Front® 2, Edge® 3, and Edge® 4
- Optional Accu-Trac Advanced Tracking Technology
- Optional WC19 crash-test approved occupied transit kit option
- Available seat sizes: 16-22" width and 14-22" depth

POWER RECLINE

- 85° - 165° of recline
- 300 lbs. weight capacity
- Available with power tilt and power adjustable seat lift
- No increase in seat-to-floor height
- Standard shear reduction and cantilever armrests

POWER POSTERIOR TILT

- 0° - 50° of tilt
- Lower seat-to-floor height
- Unmatched adjustment and durability
- Maintains center of gravity shift

DRIVE SPEEDS IN ANTERIOR TILT

- Up to 10° of anterior tilt, up to full lift: **2 MPH**
- 10-20° of anterior tilt, up to 6" of lift: **0.5 MPH**
- 10-20° of anterior tilt, over 6" of lift: **DRIVE LOCKOUT**
- Over 20° of anterior tilt: **DRIVE LOCKOUT**

COLORS

* See website for a complete list of color options.



LEMON CRUSH

RED CORVETTE

DEEP PURPLE


EVOLUTION ORANGE

RASPBERRY BERET

GREEN MACHINE

BLUE JEAN

TRU-Balance® 4
with 12" iLevel®
Power Elevation

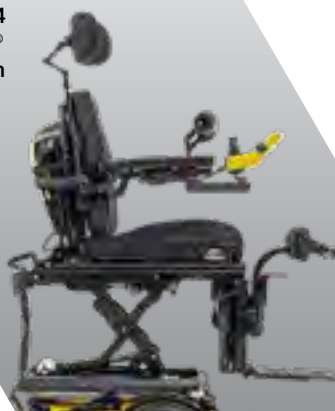

Knee Blocks



FDA Class II
Medical Device

TRU-Balance®
Flex Seating



Contoured Seating



STU LUNDY
shown in the Edge® 4
with 12" iLevel® power
seat elevation



SEATING

AND POSITIONING

TRU-BALANCE® FLEX SEATING



No-charge static seat for the Edge 3 Stretto®, Edge® 3,
Q6 Edge® 2.0, J4 and Rival®

Easily adjustable, versatile design that can grow with the user

Same ease of width and depth adjustment as TRU-Balance® 3

Seat sizes 15"-22" wide and 15"-20" deep

FDA Class II
Medical Device

CONTOURED SEATING



Black vinyl captain's seat with optional armrests

Compatible with TRU-Balance 3 power articulating foot platform

Available in high back, medium back, solid seat pan,
limited recline, and depth-adjustable



PEDIATRIC SEATING

Easily adjustable, this versatile design can grow with the user

Several width/depth adjustment packages start at 10"

Accepts a wide variety of after-market seats and backs

Available in numerous shroud colors that kids will love

Pediatric Seating




Stealth Products® CUSHIONS

FEATURING EXCLUSIVE
coolcore
CHEMICAL-FREE COOLING TECHNOLOGY



Convenient built-in side/front pouches

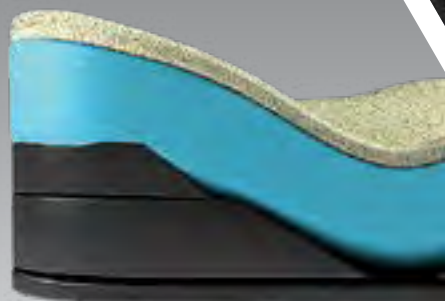

Constructed with FRG Sunmate
flame-retardant material



30+ UPF



Dry



Chemical Free



Breathable Comfort

"While the other makers of cooling fabrics use a number of chemicals to enhance wicking, transportation and evaporation, Coolcore fabrics are chemical free. To make that possible, we changed the essence of cooling fabric technology right down to its very core."

COOLCORE®



STEALTH PRODUCTS®

CUSHIONS

Stealth Products® Cushions feature patented Coolcore® technology, a chemical free, four-way stretch material that offers moisture wicking, moisture transportation, and moisture evaporation properties. Our cushions also use a color guide that identifies the suggested HCPC code and function family by a unique patch.


FEATURES





- Chemical free
- Regulates temperature
- Manages and distributes your body heat
- Ships out within 2 days on or off chair

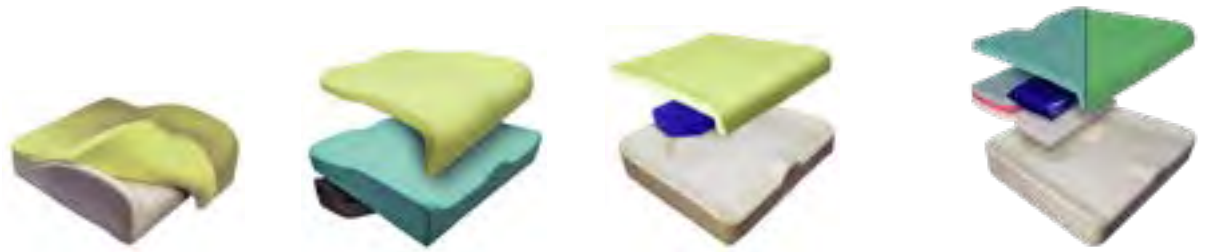

Easy-pull rubber loop zipper handles





| | TRU-Comfort 2 SPP | Essence SPP | Solution® SPP | Spectrum® Foam/Gel SPP |
|--------------------------------|---|---|---|--|
| DESCRIPTION¹ |  Multi-layered foam Deep contour, CoolCore® cushion cover has “pockets” on sides |  Low profile Tru-Comfort 2, lateral and pre-ischial shelf Contoured, HD molded foam base, coccyx and sacral relief |  Standard viscoelastic foam top with twin cell aqueous viscous gel (heavier gel) insert in ischial well, mid-level contour |  SunMate FRG foam top (SunMate FRG XXSoft (blue); Medsoft (green)) with multi-layered foam insert or twin cell ultra-lightweight viscous gel insert in ischial well, sacral relief |
| WEIGHT CAPACITY | 14-15” W: 175 lbs. 16-21” W: 350 lbs. 22+” W: 450 lbs. | 10-15” W: 175 lbs. 16-21” W: 350 lbs. 22” W: 450 lbs. | 10-15” W: 175 lbs. 16-20” W: 350 lbs. 22+” W: 500 lbs. | 14” W: 175 lbs. Foam & Gel 16-20” W: 350 lbs. Foam & Gel 22+” W: 650 lbs. Foam; 550 lbs. Gel |
| OVERALL WEIGHT | 4 lbs. | 3.5 - 5 lbs. | 5.4 lbs. | Foam: 3 lbs. Gel: 3.5 lbs. |
| HEIGHT | 4.5” | 3.5” | 3” | 3” |
| DIMENSIONS² | 14x14” up to 24x24” Odd sizes in 15-21” W SOF available | 10x10” up to 22x24” Odd sizes in 15-21” W SOF available | 10x10” up to 24x24” Odd sizes in 15-19” W SOF available | Foam: 14x18” up to 32x28” Gel: 14x18” up to 24x26” Odd sizes in 17” and 19” W SOF available |
| MSRP | \$420 - \$525 | \$420 - \$625 | \$394 - \$452 | \$483 - \$557 (\$620 for Spectrum Foam SPP larger than 24”) |
| HCPCS | E2607/E2608 | E2607/E2608 | E2607/E2608 | E2607/E2608 |
| ALLOWANCES³ | \$332.23 / \$398.98 | \$332.23 / \$398.98 | \$332.23 / \$398.98 | \$332.23 / \$398.98 |

- 1) CoolCore® is standard on all cushions except the Simplicity G.
- 2) Special order sizes available for all cushions except the Glacial SP and Simplicity G. Contact sales for additional assistance.
- 3) Allowances listed are estimations for cushions billed in conjunction with Group 3 power chair bases using the KU modifier and may vary by state.







Stealth® ADI Q-Back 2.0

FEATURES

- 450-lb. weight capacity
- CoolCore® Cover with 1.25” of Medium Soft Sunmate® FRG foam
- Two styles: AL Style (Aluminum) and AFT Style (Aluminum Flat Top)
- Available configurations: Active Contour (2.5” depth), Deep Contour (4” depth), and Extra Deep Contour (6” depth)
- Widths: 14”-22”
- Shell heights: Standard (16”), Tall (18”), Extra Tall (20”), Extra Extra Tall (22”)



| Premiere P | Glacial SP | Zen SP | Simplicity G |
|---|--|--|---|
|  <p>Gray X-Soft SunMate FRG foam topper, high density positioning wedges built in for offloading.</p> |  <p>New SIM G base, fluid resistant overlay, top and sides TPU heat-sealed, low profile, lightweight cooling gel topper</p> |  <p>New SIM G base, SunMate FRG XXSoft (blue topper), low profile, lightweight</p> |  <p>High-density foam Stealth-Tek cover standard More contour to sacral/coccyx relief well</p> |
| 14-15" W: 175 lbs. 16-21" W: 350 lbs. 22" W: 500 lbs. | 10-15" W: 175 lbs. 16-21" W: 350 lbs. 22+" W: 400 lbs. | 10-15" W: 175 lbs. 16-21" W: 350 lbs. 22+": 400 lbs. | 10-15" W: 175 lbs. 16-21" W: 350 lbs. 22+" W: 400 lbs. |
| 3.5 lbs. | 3.5 lbs. | 2.5 lbs. | 2 lbs. |
| 3.5" | 2.5" | 3" | 2" |
| 14x13" up to 22x22" Odd sizes up to 21" W SOF available | 10x10" up to 24x24" Odd sizes in 15-21" W No SOF | 10x10" up to 24x24" Odd sizes in 15-21" W SOF available | 10x10" up to 24x24" Odd sizes in 15-21" W No SOF |
| \$341 - \$499 | \$347 - \$462 | \$331 - \$446 | \$116 - \$173 |
| E2605/E2606 | E2603/E2604 | E2603/E2604 | E2601/E2602 |
| \$308.53 / \$481.33 | \$173.74 / \$215.97 | \$173.74 / \$215.97 | \$70.10 / \$136.85 |

DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement. HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. Payer coding, coverage, and bundling guidelines may apply.





FDA Class II
Medical Device



FDA Class II
Medical Device




Q-LOGIC 3 EX
ENHANCED DISPLAY



Q-LOGIC 3 EX
JOYSTICK CONTROLLER

Q-Logic 3 iAccess Programmable Switchbox
available with push buttons or toggles

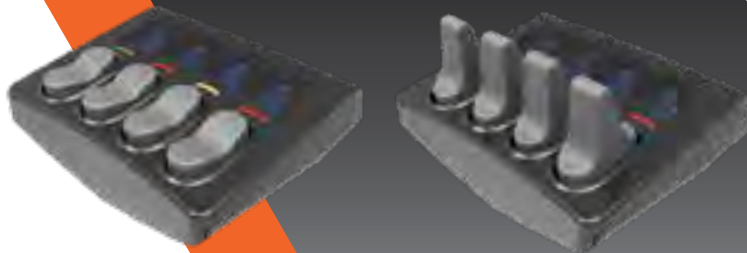

Q-Logic 3 Stand
Alone Joystick



"MY EDGE 3 GIVES ME THE OPPORTUNITY TO USE ADVANCED TECHNOLOGY THAT GIVES ME SO MUCH INDEPENDENCE."

ZOE HERNANDEZ

The Q-Logic 3e features a magnetic charging port for clients with limited dexterity



Q-LOGIC 3e
JOYSTICK
CONTROLLER

FEATURING
Bluetooth®

Q-LOGIC 3

ADVANCED DRIVE CONTROLS


The Q-Logic 3 Advanced Drive Controls System from Quantum Rehab® is like no other power wheelchair electronics system. Designed from clinician and consumer feedback, Q-Logic 3 is user-friendly and highly customizable.

FEATURES

- Wireless Bluetooth® programming through Windows devices
- Bluetooth access to computer, tablet, and phone functions
- iAccess programmable push button and toggle options for customized seating access
- Clinic Mode:
 - Quickly set up multiple input devices in client evaluations
 - Automatically assigned drive modes to input devices



Q-Logic 3
home screen

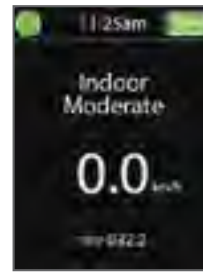

Q-Logic 3
drive screen

ENHANCED DISPLAY

The optional Q-Logic 3 EX Enhanced Display serves as a separate screen used with specialty drive controls to display system data. With the Enhanced Display's built-in digital controls, Q-Logic 3 can operate a wide range of appliances that use infrared (IR) remote controls, such as televisions, and serves as a Bluetooth® computer mouse emulator.



Q-Logic 3 Specialty Controls
Interface Module



NON-EXPANDABLE

JOYSTICK OPTIONS

SIP-N-PUFF



Q-Logic NE and NE+



Q-Logic NES



Stealth® i-Drive®
Sip & Puff

EXPANDABLE (Q-LOGIC 3)



Q-Logic 3 EX
Joystick



Q-Logic 3e
Joystick



Q-Logic 3 EX
Enhanced Display



Q-Logic 3
Stand Alone
Joystick



Q-Logic 3 iAccess
Programmable
Switchbox available
with push buttons
or toggles



Q-Logic 3
Specialty Interface
Module

PROPORTIONAL DRIVE CONTROLS



Mo-Vis™ Proportional
Head Control



Mo-Vis™
All-Round
Joystick



Mo-Vis™
All-Round Light
Joystick



Mo-Vis™
Multi Joystick



Mo-Vis™
Foot Control



Stealth® Precision
Mini Proportional
Joystick



Stealth®
Mushroom
Joystick



Mo-Vis™
Micro Joystick



Stealth® i-Drive®
Control Harness

MULTI-SWITCH



Stealth® i-Drive® Pro Series
3-Switch Head Array



Stealth® i-Drive® Tri-Array,
3-Switch Head Array



Stealth® i-Drive®
Curved Array

SINGLE SWITCH



Stealth®
Egg Switch



Mo-Vis™
Twister Basic
Button



Mo-Vis™
Twister Satellite
Button



Stealth® i-Drive®
Eclipse Tray with
Four Capacitive
Proximity Sensors



Stealth® i-Drive®
Eclipse Tray with
Four Fiber Optic
Proximity Sensors

TRAY ARRAYS

MISC. ACCESSORIES



Front and Rear
LED Fender Lights



Hydration System



Transfer Bars (Pair)



Backup Camera



Personal Item Hook



Phone Holder



Personal Items Pouch



Glove Box

NEW COLORS



INFERNO RED



ROYALE BLUE



BERMUDA BLUE



COSMIC PURPLE



FOREST GREEN



GHOST GRAY



ICE PINK



STEEL BLUE



BLACK CAMO

- COMING SOON -



NOT ALL COLOR OPTIONS ARE AVAILABLE ON ALL PRODUCTS. TO CONFIRM AVAILABILITY, VISIT [QUANTUMREHAB.COM](https://www.quantumrehab.com) OR SCAN OUR QR CODE.

LEGACY COLORS



ELECTRIC BLUE



BLUE JEAN (MATTE)



MIDNIGHT BLUE



DEEP PURPLE



RED CORVETTE



RASPBERRY BERET



PRETTY IN PINK



LEMON CRUSH



EVOLUTION ORANGE



GREEN MACHINE



GREY STREET



BACK IN BLACK



NIGHT TRAIN



MILITARY GREEN




MILITARY DIGITAL CAMO



FADED GLORY

QUANTUM® OUTBACK COLOR OPTIONS

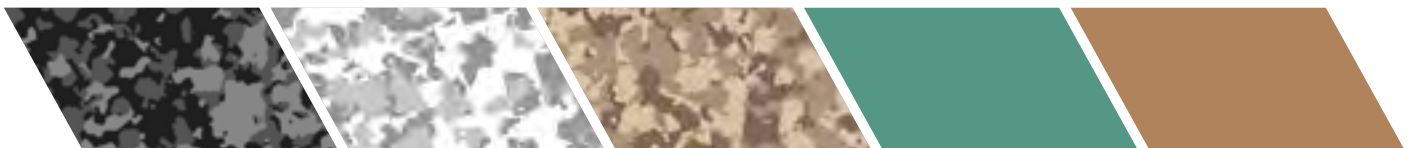

BLACK CAMO

WHITE CAMO

DESERT CAMO

CADET

TAN



INTERACTIVE ASSIST

Reduce service visits with this innovative app.

- SCREEN MIRRORING
- REAL-TIME MONITORING
- EFFICIENCY IMPROVEMENTS
- TROUBLESHOOTING
- BATTERY MANAGEMENT
- DIAGNOSTIC HISTORY
- BASIC PROGRAMMING
- REMOTE ACCESS
- CHANGE LOG HISTORY
- SAVING XPA FILES
- MOTOR USAGE HISTORY
- TEXT CHATTING
- PROGRAMMING IR MENU
- ACCESS TO BLUETOOTH FEATURES
- DEALER-TO-DEALER PROGRAMMING ACCESS

Q-Logic 3 now features Interactive Assist. Interactive Assist from Quantum Rehab® is an app which provides a direct real-time connection from a power chair's electronics to a provider's programming station. The connection provides the technician with full system and diagnostics information along with tools, such as a real-time mirror image of the power chair's electronics display, to eliminate wasteful trips.



Download Interactive Assist today!


ECON-I

Convenience at your fingertips!

The Econ-I app is now available for download on Apple® iPhone® devices, allowing you to program Quantum® Power Chairs conveniently from your iPhone or iPad. Compatible devices include the iPhone 6S/6s Plus/SE or higher, and the iPad Air/iPad Mini 2/iPad Pro/iPad 5th generation or higher. All devices must support iOS version 12 or higher.



Download ECON-I today!

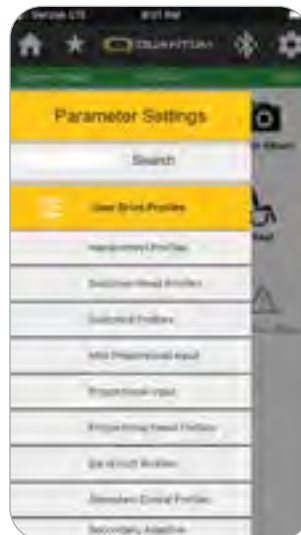


APPS

SOLUTIONS AT YOUR FINGERTIPS

PQ-365

All Quantum® and Pride® product information combined into one app!

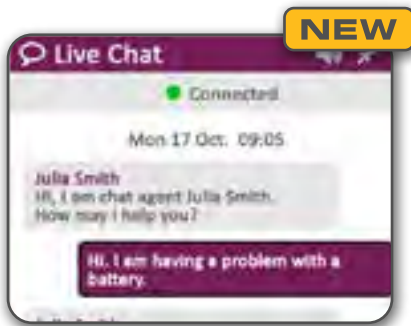

LIVE CHAT



VOICE SEARCH



ZOOM FEATURE FOR
PRODUCT IMAGES



SERIAL NUMBER AND
ORDER SEARCH



ILLUSTRATED
PARTS BREAKDOWNS



IMAGE SEND OPTION



Download PQ-365 today!

SALES SUPPORT 866.800.2002

Quantum Sales Support assists with design, ordering, pricing and quality assurance.

CUSTOM GROUP 877.465.2999

The Quantum Custom Group provides innovative custom solutions and specialized product modifications.

TECHNICAL SUPPORT EXT. 5200

Quantum Technical Support offers assistance with troubleshooting and parts order placement.

PRODUCT PLANNING & REIMBURSEMENT

Assistance with understanding Medicare codes and coverage criteria, product planning, provision guidelines, claims coding, pre-screens and assistance in regard to high-end rehab clients.

SET TO SPEC SERVICE

Set to Spec Service from Quantum Rehab® is designed to assist you in delivering Q6 Series Power Chairs to your clients. We will pre-set, to your specifications, seat-to-floor heights, seating center of gravity, armrest heights, front riggings, thoracic and/or thigh laterals.



QUANTUMREHAB.COM

Quantum's official website is the ultimate source for everything you need to know about Quantum® products. Highlights include detailed product specifications, order forms, videos, user testimonials and much more.


SERVICE

AND SUPPORT



LIVE LENS

Connect with technical service for remote, hands-on troubleshooting from anywhere in the field. Contact our Tech Service Dept. for details.

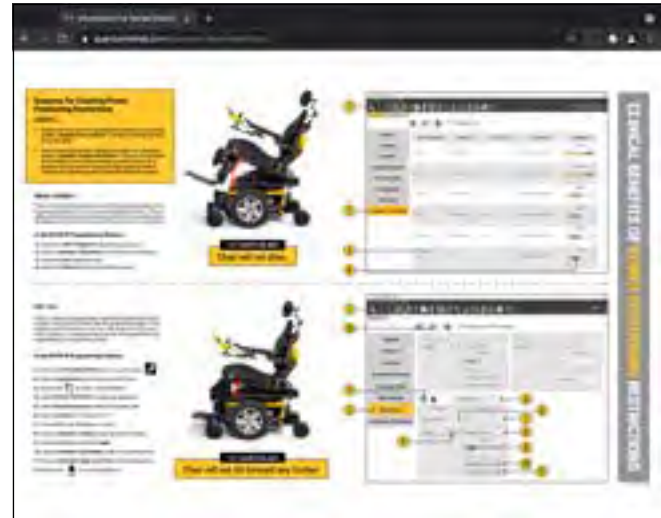

CLINICAL RESOURCES

The Clinical Resources section of the Quantum website is specifically tailored to clinicians' needs, featuring education, order forms, clinical references, and reimbursement information.



SOCIAL MEDIA

Stay up to date on the latest Quantum® events and product releases through Twitter, Facebook, YouTube, Instagram and TikTok.



STEALTHPRODUCTS.COM

Stealth Products® launched in 1999, in response to a need for versatile, high-quality head and neck positioning systems, as well as switch site location hardware. Their products include: head and body support/positioning, joystick mounts, trays, seating, cushions, and pediatric mobility products.



QUANTUM®
#1 FOR REHAB POWER

- 1) Standard
- 2) Optional
- 3) Requires HD Suspension
- 4) Speed varies with user weight, terrain type, battery charge, battery condition and tire pressure.
- 5) Without front riggings
- 6) Battery capacity may change when adding power positioning systems.
- 7) Weight capacity is dictated by the selection of seating.
- 8) All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 300 lbs. weight capacity and solid tires (add 3/8" with pneumatic tires). Depths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.
- 9) 250 lbs. weight capacity with TRU-Balance®
- 10) A minimum seating size of 20" x 20" is required with iLevel® on the Q6 Edge® HD; maximum seating size is 24" x 24"
- 11) Without foot riggings
- 12) All Quantum units can be ordered with specialty controls
- 13) 6.25 mph speed available only with Q-Logic 3.
- 14) Battery weight may vary +/- 3 lbs. based on manufacturer.
- 15) Reduced to 400 lbs. with iLevel® power adjustable seat height
- 16) Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 0.25" lower when the unit is loaded.
- 17) 7 mph can only be achieved with optional 4-pole motors and Q-Logic 3. Maximum speed without 4-pole motors is 6.25 mph with Q-Logic 3.

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement. The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.

Quantum® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.



| | EDGE® 4 | EDGE 3 STRETTO® |
|--|--|--|
| HCPCS CODES For more information on HCPCS code coverage, please consult our Product Planning & Reimbursement Center. | Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861) | Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861) |
| AVAILABLE ON VA CONTRACT | Coming Soon | FSS Contract V797D-70001 |
| WEIGHT CAPACITY ⁷ | 300 lbs. | 300 lbs. |
| MAXIMUM SPEED ⁴ | Up to 7 mph (4.5 mph with iLevel®) | Up to 6.25 mph (3.5 mph with iLevel®) ¹³ |
| GROUND CLEARANCE | 2.87" | 2.9" |
| TURNING RADIUS ⁵ | 22.8" | 21" |
| LENGTH ¹¹ | 38.27" | 36" |
| WIDTH | 24.02" | 21.5" (14" drive wheels) 20.47" (12.5" drive wheels) |
| DRIVE WHEELS | 14" pneumatic or flat free tires Optional: 14" knobby | 14" or 12.5" pneumatic or flat free tires; Optional: 14" knobby |
| CASTER WHEELS | 7" flat free | Front: 5" flat free Rear: 6" flat free |
| ANTI-TIP WHEELS | N/A | N/A |
| BASE WEIGHT | 154 lbs. without batteries | 126.4 lbs. |
| BATTERY WEIGHT ¹⁴ | 47.8 lbs. each | 40 AH: 31 lbs. each; NF-22: 40 lbs. each |
| DRIVE WHEEL CONFIGURATION | Mid-wheel Drive | Mid-wheel Drive |
| MOTOR PACKAGES | 4-pole | 2-pole |
| BATTERIES ⁶ | Group 24 (C5=65AH; C20=75AH) | 40AH (2) or NF-22 (2) |
| BATTERY CHARGER | 8A Off-board | 8A Off-board |
| AVAILABLE ELECTRONICS ¹² | 120A Q-Logic 3 90A Q-Logic NES | 120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+ |
| SEATING SIZES ⁸ | | |
| TRU-BALANCE® 4 ANTERIOR TILT | W: 16-22" D: 14-22" | W: 16-22" D: 14-22" |
| TRU-BALANCE® 4 POWER TILT | W: 12-22" D: 12-22" | W: 12-22" D: 12-22" |
| SYNERGY® | N/A | W: 10-22" D: 10-22" |
| TRU-BALANCE® 3 STATIC | W: 12-22" D: 12-20" | W: 12-22" D: 12-22" |
| TRU-BALANCE® 3 POWER TILT | W: 12-22" D: 12-22" | W: 12-22" D: 12-22" |
| TRU-BALANCE® 3 LIFT & TILT | W: 12-22" D: 12-22" | W: 12-22" D: 12-22" |
| TRU-BALANCE® 3 HD STATIC | N/A | N/A |
| TRU-BALANCE® 3 HD POWER TILT | N/A | N/A |
| TRU-BALANCE® 3 HD LIFT & TILT | N/A | N/A |
| TRU-BALANCE® 2 HD STATIC | N/A | N/A |
| TRU-BALANCE® 2 HD POWER TILT | N/A | N/A |
| TRU-BALANCE® FLEX SEAT | W: 10-22" D: 10-20" | W: 15-22" D: 15-20" |
| Q4 POWER TILT | N/A | N/A |
| SEAT-TO-FLOOR HEIGHTS ^{8,13,16} | | |
| TRU-BALANCE® 4 ANTERIOR TILT | 18.25-21.25" | 18-19" |
| TRU-BALANCE® 4 POWER TILT | 17.75-19.75" | 17.5-19.5" |
| SYNERGY® / STATIC | N/A | 16"-18" (Static only) |
| TRU-BALANCE® 3 POWER TILT | 17.75-19.75" | 17.5-19.5" |
| LIFT & TILT | 17.5-19.5" | 18-20" |
| TRU-BALANCE® FLEX SEAT | 16.25-18.25" | 16.5-18.5" |
| Q4 SEATING | N/A | N/A |
| Q CAPTAIN SEAT | SSP: 16.5-18.5" (Flex); 18-20" (Lift) Captain: 20.5-22.5" (Flex); 22-24" (Lift) | SSP: 16.5-18.5" (Flex); 18-20" (Lift) Captain: 20.5-22.5" (Flex); 22-24" (Lift) |
| OPTIONAL POWER ELEVATION | 12" iLevel® at 4.5 mph | 12" iLevel® at 3.5 mph |



| EDGE® 3 | | EDGE® 2.0 | | EDGE® HD | | OUTBACK | |
|--|--|--|--|--|--|---|--|
| Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861) | | Group 3: Standard (K0848/K0849), Single Power (K0856/K0857), Multiple Power (K0861) | | Group 3: Heavy Duty (K0850/K0851), Single Power HD (K0858), Multiple Power HD (K0862) | | K0899 - Power Mobility Device, Not coded by DME PDAC | |
| FSS Contract V797D-70001 | | FSS Contract V797D-70001 | | FSS Contract V797D-70001 | | Coming Soon to FSS Contract | |
| 300 lbs. | | 300 lbs. | | 450 lbs. ¹⁵ | | 300 lbs. (350 lbs. w/ Q-Captain) | |
| Up to 6.25 mph ¹³ Optional: Up to 7 mph ¹⁷ (4.5 mph with iLevel®) | | Up to 6.25 mph ¹³ (3.5 mph with iLevel®) | | Up to 4.7 mph (3.5 mph with iLevel®) | | Up to 5 mph (6 mph with high speed package) | |
| 2.7" (battery tray) | | 2.7" (battery tray) | | 3" (base frame) | | 4.92" (at motors) | |
| 20.5" | | 20.5" | | 19.75" | | 49.2" | |
| 35.5" | | 35.5" | | 35.6" | | 42" | |
| 24.4" | | 24.4" | | 26.25" | | 28.82" | |
| 14" Solid Optional: 14" Pneumatic; 14" Knobby | | 14" Solid Optional: 14" Pneumatic; 14" Knobby | | 14" Run Flat Optional: 14" Pneumatic; 14" Knobby | | 15" pneumatic | |
| Front: 6" Solid Rear: 6" Solid | | Front: 6" Solid Rear: 6" Solid | | Front: 6" Solid Rear: 6" Solid | | N/A | |
| N/A | | N/A | | N/A | | N/A | |
| 138.6 lbs. | | 138.6 lbs. | | 143.5 lbs. | | 168.6 lbs. | |
| 40 lbs. each | | 40 lbs. each | | 52 lbs. each | | 58.6 lbs. | |
| Mid-wheel Drive | | Mid-wheel Drive | | Mid-wheel Drive | | 4 Motor, Independent Drive | |
| 2-pole (Optional: 4-pole) | | 2-pole | | 4-pole | | 4-pole | |
| NF-22 (2) | | NF-22 (2) | | Group 24 (2) | | Group 27 | |
| 8A Off-board | | 8A Off-board | | 8A Off-board | | 8A Off-board | |
| 120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+ | | 120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+ | | 120A Q-Logic 3; 75A Q-Logic NE; 75A Q-Logic NE+ | | 120A Q-Logic 3 | |
| W: 16-22" D: 14-22" | | N/A | | N/A | | N/A | |
| W: 12-22" D: 12-22" | | N/A | | N/A | | N/A | |
| W: 10-22" D: 10-22" | | W: 10-22" D: 10-22" | | W: 20-28" D: 15-28" | | N/A | |
| W: 12-22" D: 12-22" | | W: 12-22" D: 12-22" | | N/A | | N/A | |
| W: 12-22" D: 12-22" | | W: 12-22" D: 12-22" | | N/A | | W: 16-22" D: 16-20" | |
| W: 12-22" D: 12-22" | | W: 12-22" D: 12-22" | | N/A | | W: 16-22" D: 16-20" | |
| N/A | | N/A | | W: 20-32" D: 20-25" | | N/A | |
| N/A | | N/A | | W: 20-32" D: 20-25" | | N/A | |
| N/A | | N/A | | W: 20-24" D: 20-25" ¹⁰ | | N/A | |
| N/A | | N/A | | N/A | | N/A | |
| N/A | | N/A | | N/A | | N/A | |
| W: 15-22" D: 15-20" | | W: 15-22" D: 15-20" | | N/A | | N/A | |
| N/A | | N/A | | N/A | | W: 16-22" D: 16-22" | |
| 17.5-19.25" | | N/A | | N/A | | N/A | |
| 17.5-19.5" | | N/A | | N/A | | N/A | |
| Synergy®: 15.875 - 20.875" Static: 16.5 - 18.5" | | Synergy®: 15.875 - 20.875" Static: 16.5 - 18.5" | | 16.5-18.5" | | N/A | |
| 17.5-19.5" | | 17.5-19.5" | | 18.5-20.5" | | N/A | |
| 17.5-19.5" | | 17.5-19.5" | | 18.5-20.5" | | N/A | |
| 16.5-18.25" | | 16.5-18.25" | | N/A | | N/A | |
| N/A | | N/A | | N/A | | N/A | |
| SSP: 16.5-18.25" (Flex); 17.5-19.5" (Lift) Captain: 20.5-22.5" (Flex); 21.5-23.25" (Lift) | | SSP: 16.5-18.25" (Flex); 17.5-19.5" (Lift) Captain: 20.5-22.5" (Flex); 21.5-23.25" (Lift) | | SSP: 16.5-18.5" (Static TB3) Captain: 20.5-22.5" (Static TB3) | | Captain: 19.5-21.5" | |
| 12" iLevel® at 4.5 mph | | 12" iLevel® at 3.5 mph | | 12" iLevel® at 3.5 mph | | 10" Safe Seat Elevation at 3.2 mph | |



QUANTUM®
#1 FOR REHAB POWER

- 1) Standard
- 2) Optional
- 3) Requires HD Suspension
- 4) Speed varies with user weight, terrain type battery charge, battery condition and tire pressure.
- 5) Without front riggings
- 6) Battery capacity may change when adding power positioning systems.
- 7) Weight capacity is dictated by the selection of seating.
- 8) All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 300 lbs. weight capacity and solid tires (add 3/8" with pneumatic tires). Depths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.
- 9) 250 lbs. weight capacity with TRU-Balance®
- 10) A minimum seating size of 20" x 20" is required with iLevel® on the Q6 Edge® HD; maximum seating size is 24" x 24"
- 11) Without foot riggings
- 12) All Quantum units can be ordered with specialty controls
- 13) 6.25 mph speed available only with Q-Logic 3.
- 14) Battery weight may vary +/- 3 lbs. based on manufacturer.
- 15) Reduced to 400 lbs. with iLevel® power adjustable seat height
- 16) Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 0.25" lower when the unit is loaded.
- 17) 7 mph can only be achieved with optional 4-pole motors and Q-Logic 3. Maximum speed without 4-pole motors is 6.25 mph with Q-Logic 3.

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement. The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.

Quantum® Power Chairs are FDA Class II Medical Devices designed to aid individuals with mobility impairments.



| | J4E | J4E HD |
|--|--|------------------------------------|
| HCPCS CODES For more information on HCPCS code coverage, please consult our Product Planning & Reimbursement Center. | Group 2: Single Power (K0835) | Group 2: Single Power (K0837) |
| AVAILABLE ON VA CONTRACT | Coming Soon | Coming Soon |
| WEIGHT CAPACITY ⁷ | 300 lbs. | 400 lbs. |
| MAXIMUM SPEED ⁴ | Up to 5 mph | Up to 5 mph |
| GROUND CLEARANCE | 3" (battery tray) | 2.7" (battery tray) |
| TURNING RADIUS ⁵ | 20.5" | 23" |
| LENGTH ¹¹ | 35.5" | 35.5" |
| WIDTH | 24.4" | 24.4" |
| DRIVE WHEELS | 14" Solid | 14" Solid |
| CASTER WHEELS | Front: 6" Solid Rear: 6" Solid | Front: 6" Solid Rear: 6" Solid |
| ANTI-TIP WHEELS | N/A | N/A |
| BASE WEIGHT | 138.6 lbs. | 138.6 lbs. |
| BATTERY WEIGHT ¹⁴ | 38 lbs. each | 38 lbs. |
| DRIVE WHEEL CONFIGURATION | Mid-wheel Drive | Mid-wheel Drive |
| MOTOR PACKAGES | 2-pole | 4-pole |
| BATTERIES ⁶ | NF-22 (2) | NF-22 |
| BATTERY CHARGER | 8A Off-board | 8A Off-board |
| AVAILABLE ELECTRONICS ¹² | 120A Q-Logic 3E 75A Q-Logic NE+ (coming soon) | 120A Q-Logic 3e 75A Q-Logic NES |
| SEATING SIZES | | |
| TRU-BALANCE® ⁴ | N/A | N/A |
| SYNERGY® | N/A | N/A |
| TRU-BALANCE® 3 STATIC | N/A | N/A |
| TRU-BALANCE® 3 POWER TILT | N/A | N/A |
| TRU-BALANCE® 3 LIFT & TILT | N/A | N/A |
| TRU-BALANCE® 3 HD STATIC | N/A | N/A |
| TRU-BALANCE® 3 HD POWER TILT | N/A | N/A |
| TRU-BALANCE® 3 HD LIFT & TILT | N/A | N/A |
| TRU-BALANCE® 2 HD STATIC | N/A | N/A |
| TRU-BALANCE® 2 HD POWER TILT | N/A | N/A |
| TRU-BALANCE® FLEX SEAT | N/A | N/A |
| Q4 POWER TILT | W: 16-20" D: 16-22" | W: 20-24" D: 20-24" |
| SEAT-TO-FLOOR HEIGHTS^{13,16} | | |
| TRU-BALANCE® ⁴ | N/A | N/A |
| SYNERGY® / STATIC ⁸ | N/A | N/A |
| POWER TILT ⁸ | N/A | N/A |
| LIFT & TILT ⁸ | N/A | N/A |
| TRU-BALANCE® FLEX SEAT ⁸ | N/A | N/A |
| Q4 POWER TILT | 17.5-19.5" | 17.5-19.5" |
| Q4 LIFT AND TILT | 17.5-19.5" | 17.5-19.5" |
| Q CAPTAIN SEAT | N/A | N/A |
| OPTIONAL POWER ELEVATION | 10" Safe Seat Elevation | 10" Safe Seat Elevation |

quantumrehab.com •     
The information contained herein is correct at the time of publication.
We reserve the right to alter specifications without prior notice.

life.quantumrehab.com
(US) 866-800-2002 • (Canada) 888-570-1113



1 F O R R E H A B P O W E R

Quantum® Power Chairs are FDA Class II Medical Devices
designed to aid individuals with mobility impairments.