

# Quantisal™ II

ORAL FLUID COLLECTION DEVICE

**The best of Quantisal in a split specimen device.**

Offering you the simplicity of collecting **two samples with a single device** that meets SAMSHA<sup>1</sup> and DOT<sup>2</sup> guidelines.



For *in vitro* diagnostic use.

[www.redplanettesting.com](http://www.redplanettesting.com)  
610-509-0567  
[info@redplanettesting.com](mailto:info@redplanettesting.com)

# Quantisal™ II

ORAL FLUID COLLECTION DEVICE

## Making split specimen collections simple, fast and convenient

While there are a number of matrices for drug testing, oral fluid testing provides a better recent-use indicator than urine<sup>3</sup> and protects against specimen adulteration and tampering—allowing you to make informed decisions to maintain a drug-free workplace.

Quantisal II is the first commercial oral fluid device with split specimen capabilities—continuing Abbott’s history of innovative solutions. Test with confidence with Quantisal.



### ORAL FLUID BENEFITS

- Non-invasive, gender-neutral specimen collections
- Collect samples anywhere
- Observed collection protects from specimen tampering or adulteration
- Excellent indication of recent drug use<sup>3</sup>
- Parent drugs present in higher concentration than metabolites<sup>4</sup>
- Easier implementation of random testing




### WHY CHOOSE QUANTISAL II?

- Meets SAMHSA Mandatory Guidelines for Federal Workplace Drug Testing Programs – Oral Fluid and DOT rule CFR 49, Part 40
- Allows simultaneous collection of A and B specimens
- Provides two 4 mL samples containing 1 mL of oral fluid and 3 mL of preservative buffer – giving adequate volume for screening and confirmation testing
- Proprietary FDA 510(k) cleared device

# How It Works


1

Collect the oral fluid sample

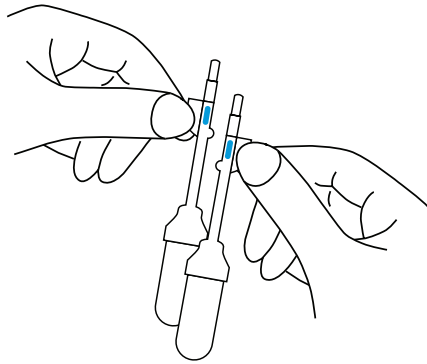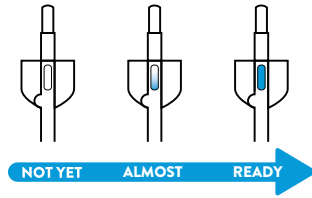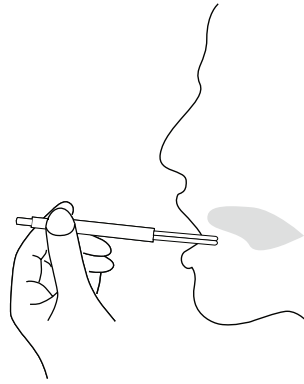

2

Hold the device by its “wings” and pull apart into two split pads/stems

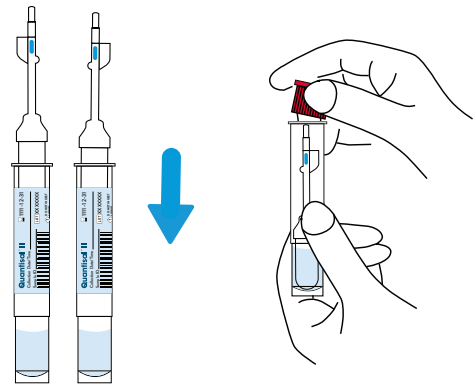

3

Insert separated pad/stem into tubes, snap on the caps and prepare shipping




The Quantisal™ II Oral Fluid Collection Device is intended for the collection, preservation and transportation of oral fluid specimens for drug of abuse testing. This device is for prescription use only.

Quantisal II device performance has been supported using the representative drug analytes in the table below:

6-AM (heroin metabolite)	Methamphetamine
Amphetamine	Nordiazepam
Buprenorphine	Opiates (including codeine and morphine)
Cocaine / Benzoylcegonine	Oxycodone
Hydrocodone	Phencyclidine (PCP)
Marijuana (THC)	Tramadol
Methadone	

The **advantages** of oral fluid.  
The **confidence** of **Quantisal**.

1. Substance Abuse and Mental Health Services Administration (SAMHSA)
2. Department of Transportation (DOT), Federal Aviation Administration, 14 CFR Part 120, Procedures for Transportation Workplace Drug and Alcohol Testing Programs: Addition of Oral Fluid Specimen Testing for Drugs, Federal Register Vol. 88 No. 84, May 02, 2023
3. White RM, Moore CM. *Detection of Drugs and Their Metabolites in Oral Fluid*. Elsevier; 2018.
4. Dept of Health and Human Services (HHS) and Substance Abuse and Mental Health Services Administration (SAMHSA), Clinical Drug Testing in Primary Care. Technical Assistance Publication Series 32. HHS Publication 12-4668, 2012

 [www.redplanettesting.com](http://www.redplanettesting.com)  
 610-509-0567  
 [info@redplanettesting.com](mailto:info@redplanettesting.com)

